

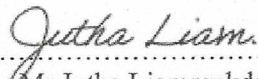
Report Test

M/C : Henan Sanqiang Medical Equipment Liability Co. Ltd
Model : SQ-H220
Serial No : HNSQ2020050063

LOT No : 00954
Date : 21 March 2025

Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
Humidification			
Humidity	50.0-90.0%	20.2 %	
Sterilization			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
Gas Injection			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
Exposure			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
Exhaust			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
Storage			
Storage Time	≥ 48 hour	≥ 48 hour	
Chemical Test			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
Biological Test			
EZ Test Biological Indicator	Green	Green	Pass


(.....)
Mr. Jutha Liam mukda
Research and Development Engineer

			LOT	Color	BOX	TOTAL / PCS
		205593				
1	Z-INDBL/SC2-1.6-10 0-5/8-CR	STERILE SURGICAL STEEL INTERNALLY THREADED DOUBLE JEWEL NAVEL BANANA	LOT NO 00954			5
		207061				
2	Z-TIMBN-1.2-10 0-3-HP	STERILE TITANIUM INTERNALLY THREADED CURVED MICRO BARBELL	LOT NO 00954			10
		207332				
1	Z-TBL-1.6-36 0-4-HP	STERILE TITANIUM BARBELL	LOT NO 00954			6
2	Z-TMBN-1.2-12 0-2.5-HP	STERILE TITANIUM CURVED MICRO BARBELL				15
		207333				
1	Z-TJNO-1.0 6.5 2.35 HP CR	STERILE TITANIUM PREMIUM CRYSTAL JEWELLED CURVED NOSE STUD	LOT NO 00954			4
2	Z-TMCB-1.2-10 0-3-HP	STERILE TITANIUM MICRO CIRCULAR BARBELL	LOT NO 00954			4
3	Z-TMLB-1.2-8 0-3-HP	STERILE TITANIUM MICRO LABRET	LOT NO 00954			8
		207341				
3	Z-TMLB-PIN-1.2-7 0-4B-HP	STERILE TITANIUM MICRO LABRET PIN	LOT NO 00954			50
		207343				
1	Z-BNDJ-1.6 12 0 5/8 CR	STERILE SURGICAL STEEL DOUBLE JEWELLED NAVEL BANANA	LOT NO 00954			160
2	Z-MLB-1.2-3 0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 00954			160
		207344				
3	Z-TMLB-1.2-10 0-3/4B-HP	STERILE TITANIUM MICRO LABRET	LOT NO 00954			50
		207357				
5	Z-BKMCB-1.2-8 0-3	STERILE BLACK PVD COATED SURGICAL STEEL MICRO CIRCULAR BARBELL	LOT NO 00954			4
6	Z-BKMLB-1.2-6 0-3	STERILE BLACK PVD COATED SURGICAL STEEL MICRO LABRET	LOT NO 00954			4
7	Z-BKMLBJ-1.2 6 0 2.5 CR	STERILE BLACK PVD COATED SURGICAL STEEL JEWELLED MICRO LABRET	LOT NO 00954			4
8	Z-DL-1.6-10 0-5	STERILE SURGICAL STEEL BARDELL	LOT NO 00954			6
9	Z-BL-1.6-20 0-5	STERILE SURGICAL STEEL BARBELL	LOT NO 00954			6
10	Z-BL-1.6-22 0-5	STERILE SURGICAL STEEL BARBELL	LOT NO 00954			4
11	Z-GPMBN-1.2-8 0-3	STERILE GOLD PVD COATED SURGICAL STEEL CURVED MICRO BARBELL	LOT NO 00954			5
12	Z-MBN-1.2-10 0-3	STERILE SURGICAL STEEL CURVED MICRO BARBELL	LOT NO 00954			4
13	Z-MBN-1.2-6 0-3	STERILE SURGICAL STEEL CURVED MICRO BARBELL	LOT NO 00954			4
14	Z-MBN-1.2-8 0-3	STERILE SURGICAL STEEL CURVED MICRO BARBELL	LOT NO 00954			4
15	Z-MLB-1.2-10 0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 00954			10
16	Z-MLBJ-1.2-10 0-3-CR	STERILE SURGICAL STEEL JEWELLED MICRO LABRET	LOT NO 00954			8
17	Z-RGMCB-1.2-8 0-3	STERILE ROSE GOLD PVD COATED SURGICAL STEEL MICRO CIRCULAR BARBELL	LOT NO 00954			4
18	Z-RGMLBJ-1.2-6 0-3-CR	STERILE ROSE GOLD PVD COATED SURGICAL STEEL JEWELLED MICRO LABRET	LOT NO 00954			5
19	Z-RGMLBJ-1.2-8 0-3-CR	STERILE ROSE GOLD PVD COATED SURGICAL STEEL JEWELLED MICRO LABRET	LOT NO 00954			5
20	Z-TINBL-1.6-12 0-5-HP	STERILE TITANIUM INTERNALLY THREADED BARBELL	LOT NO 00954			2
21	Z-TINBL-1.6-20 0-5-HP	STERILE TITANIUM INTERNALLY THREADED BARBELL	LOT NO 00954			2
22	Z-TINBL-1.6-22 0-5-HP	STERILE TITANIUM INTERNALLY THREADED BARBELL	LOT NO 00954			2
23	Z-TINMLB-1.2-8 0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO LABRET	LOT NO 00954			6
24	Z-TMCN-1.2-12 0-3-HP	STERILE TITANIUM MICRO CURVED BARDELL	LOT NO 00954			2
25	Z-TMCN-1.2-6 0-3-HP	STERILE TITANIUM CURVED MICRO BARBELL	LOT NO 00954			4
26	Z-TMCN-1.2-8 0-3-HP	STERILE TITANIUM CURVED MICRO BARBELL	LOT NO 00954			4
27	Z-TMCB-1.2-12 0-3-HP	STERILE TITANIUM MICRO CIRCULAR BARBELL	LOT NO 00954			4
28	Z-TMCB-1.2-8 0-3-HP	STERILE TITANIUM MICRO CIRCULAR BARBELL	LOT NO 00954			4
29	Z-ZGMLBJ-1.2-10 0-3-CR	STERILE ZIRCON GOLD PVD SURGICAL STEEL JEWELLED MICRO LABRET	LOT NO 00954			8
30	Z-TMCB-1.2-12 0-3-HP	STERILE TITANIUM MICRO CIRCULAR BARBELL	LOT NO 00954			4
		2534830				
7	PT-Z-GPTMLBJ-1.2-8 0-3-CR	STERILE GOLD PVD COATED TITANIUM PREMIUM CRYSTAL JEWELLED MICRO LABRET	LOT NO 00954			150
		2536556				
3	PT-Z-GPTMCB-1.2-10 0-3	STERILE GOLD PVD COATED TITANIUM MICRO CIRCULAR BARBELL	LOT NO 00954			50
5	PT-Z-GPTMLBJ-1.2-8 0-3-CR	STERILE GOLD PVD COATED TITANIUM PREMIUM CRYSTAL JEWELLED MICRO LABRET	LOT NO 00954			100
7	PT-Z-TBL-1.6-34 0-5	STERILIZED TITANIUM BARDELL FOR PTIT ROC	LOT NO 00954			150
10	PT-Z-TMLBJ-1.2-10 0-3-CR	STERILE TITANIUM MICRO LABRET WITH JEWELLED BALL FOR PTIT ROC	LOT NO 00954			100
18	PT-Z-TMLB-1.2-8 0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 00954		SET/150	1500
		2536567				
1	QL-Z-BKMLB-1.2-8 0-3	STERILE BLACK PVD COATED SURGICAL STEEL MICRO BARDELL FOR QUILLIAN	LOT NO 00954			100
5	QL-Z-MBL-1.2-10 0-3	STERILE SURGICAL STEEL MICRO BARDELL FOR QUILLIAN	LOT NO 00954			50
6	QL-Z-MLB-1.2-6 0-3	STERILE SURGICAL STEEL MICRO LABRET FOR QUILLIAN	LOT NO 00954			100
						2821

Sterilization:
Temperature: 44.4°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 19:19

Add EO gas
Temperature: 44.5°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 19:28

Evacuation
Temperature: 44.3°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 19:34

Keeping Temp
Temperature: 43.9°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 19:34

Free pressure qualified
Temperature: 42.7°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 19:32

Evacuation
Temperature: 41.9°C
Humidity: 20.20%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 19:39

Heat up
Temperature: 42.7°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 19:25

Setting:
Temperature: 45.0°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 19:25

21/3/25
101954

Air replacement:
Temperature: 45.6°C
Humidity: 20.20%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 1
round time: 15min(s)
2025-3-21 17:46

Air replacement:
Temperature: 45.4°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 2
round time: 15min(s)
2025-3-21 17:58

Air replacement:
Temperature: 45.4°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 3
round time: 15min(s)
2025-3-21 17:56

Air replacement:
Temperature: 45.2°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 4
round time: 15min(s)
2025-3-21 18:45

Air replacement:
Temperature: 45.4°C
Humidity: 20.20%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 18:29

Sterilization:
Temperature: 44.4°C
Humidity: 20.20%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 18:29

Sterilization:
Temperature: 44.5°C
Humidity: 20.18%
Pressure: -1.60Pa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-3-21 18:21

