

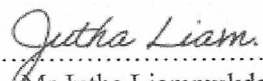
Report Test

M/C : Henan Sanqiang Medical Equipment Liability Co. Ltd
 Model : SQ-H220
 Serial No : HNSQ2020050063

LOT No : 00958
 Date : 27 March 2025

Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
Humidification			
Humidity	50.0-90.0%	20.2 %	
Sterilization			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
Gas Injection			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
Exposure			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
Exhaust			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
Storage			
Storage Time	≥ 48 hour	≥ 48 hour	
Chemical Test			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
Biological Test			
EZ Test Biological Indicator	Green	Green	Pass


 (.....)
 Mr. Jutha Liammukda
 Research and Development Engineer

			LOT	Color	BOX	TOTAL / PCS
		207627				
1	Z-XINLBTJD-1.2-12.0-3-CL-CR	STERILE BIOFLEX INTERNAL LABRET WITH JEWELLED TITANIUM DISK	LOT NO 00958			10
		207573				
2	Z-BN-1.6-16.0-4	STERILE SURGICAL STEEL CURVED BARBELL	LOT NO 00958			100
		207459				
2	Z-TENDJ-1.6-12.0-5/8-HP-CR	STERILE TITANIUM DOUBLE JEWELLED NAVEL BANANA	LOT NO 00958			50
		207434				
1	Z-TMBN-PIN-1.2-12.0-HP	STERILE TITANIUM MICRO CURVED BARBELL PIN	LOT NO 00958			30
2	Z-TMBN-PIN-1.2-7.0-HP	STERILE TITANIUM MICRO CURVED BARBELL PIN	LOT NO 00958			30
3	Z-TMBN-PIN-1.2-8.0-HP	STERILE TITANIUM MICRO CURVED BARBELL PIN	LOT NO 00958			30
4	Z-TMBN-PIN-1.2-10.0-HP	STERILE TITANIUM MICRO CURVED BARBELL PIN	LOT NO 00958			30
		207408				
13	Z-RGSCRC007-1.2-8.0-CR	STERILE ROSE GOLD PVD COATED SURGICAL STEEL PRONG SET JEWELLED ROOK CLICKER	LOT NO 00958			1
		207379				
1	Z-GPBL-1.6-14.0-5	STERILE GOLD PVD COATED SURGICAL STEEL BARBELL	LOT NO 00958			50
4	Z-MLB-1.2-9.0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 00958			150
		2534572				
1	Z-G*BTBNDJ-OP-1.6-10.0-5/8-CR	STERILE GOLD PVD 18K COATED TITANIUM DOUBLE OPTIMA CRYSTAL JEWELLED NAVEL BANANA	LOT NO 00958			20
		2534630				
6	PT-Z-G*TMLB-1.2-8.0-3	STERILE GOLD PVD COATED TITANIUM MICRO LABRET	LOT NO 00958			500
		2535363				
10	PT-Z-TBN-1.6-10.0-4	STERILE TITANIUM TITANIUM CURVED BARBELL FOR PTIT ROC	LOT NO 00958			50
		2535615				
1	TDI-Z-MLB-OP-1.2-10.0-3-CR	STERILE SURGICAL STEEL CZECH CRYSTAL JEWELLED MICRO LABRET FOR TDI	LOT NO 00958			300
		2536566				
4	PT-Z-G*TMLB-1.2-8.0-3	STERILE GOLD PVD COATED TITANIUM MICRO LABRET	LOT NO 00958			400
9	PT-Z-TLCHC-0.6-5.2-HP	STERILE TITANIUM LARGE CURVED BALL NOSE STUD	LOT NO 00958			100
16	PT-Z-TL*EL-1.6-16.0-5	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 00958		SET/25	250
		2537579				
1	PT-Z-TIMMLB-1.2-6.0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO LABRET FOR PTIT ROC	LOT NO 00958			10
		2537635				
1	Z-LJNC-S-0.6-19.0-2.35-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELLED STRAIGHT NOSE STUD 19MM	LOT NO 00958			20
		2537667				
1	PT-Z-GPT*CSH-1.2-10.0	STERILE GOLD PVD COATED TITANIUM HINGED SEGMENT RING	LOT NO 00958			200
		2536905				
3	SV-STERILIZE	SERVICE STERILE	LOT NO 00958			10
		2536792				
3	SV-STERILIZE	SERVICE STERILE				10
						2301

28/9/25 104958 EMI Air replacement: Temperature: 45.5t Humidity: 20.2RH Pressures: -1.0kPa Sterilization time: 3hour Air replace rounds: 1 round time: 15min(s) 2025-9-27 17:39 Air replacement: Temperature: 45.1t Humidity: 20.1RH Pressures: -1.0kPa Sterilization time: 3hour Air replace rounds: 2 round time: 15min(s) 2025-9-27 17:19 Air replacement: Temperature: 45.5t Humidity: 20.2RH Pressures: -1.0kPa Sterilization time: 3hour Air replace rounds: 3 round time: 15min(s) 2025-9-27 16:50 Air replacement: Temperature: 45.3t Humidity: 20.2RH Pressures: -1.0kPa Sterilization time: 3hour Air replace rounds: 4 round time: 15min(s) 2025-9-27 16:30 Air replacement: Temperature: 45.1t Humidity: 20.2RH Pressures: -1.0kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 16:19 Sterilization: Temperature: 44.7t Humidity: 20.2RH Pressures: -0.8kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 16:12 Sterilization: Temperature: 44.6t Humidity: 20.2RH Pressures: -0.9kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 13:14	Sterilization: Temperature: 44.8t Humidity: 20.1RH Pressures: -0.9kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 13:12 Add EO gas Temperature: 44.5t Humidity: 20.1RH Pressures: -4.0kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 13:12 Evacuation: Temperature: 44.4t Humidity: 20.1RH Pressures: -0.8kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 13:06 Keeping Temp Temperature: 44.2t Humidity: 20.1RH Pressures: -1.0kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 12:06 Keep pressure qualified Temperature: 43.0t Humidity: 20.1RH Pressures: -55.0kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 12:04 Evacuation: Temperature: 42.0t Humidity: 20.1RH Pressures: -0.8kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 11:53 Heating up: Temperature: 25.1t Humidity: 20.1RH Pressures: -0.8kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 11:26 Setting: Temperature: 45.0t Humidity: 55.0RH Pressures: -50.0kPa Sterilization time: 3hour Air replace rounds: 5 round time: 15min(s) 2025-9-27 11:20
--	---