

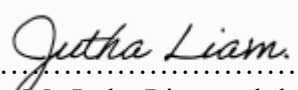
Report Test

M/C : Henan Sanqiang Medical Equipment Liability Co. Ltd
Model : SQ-H220
Serial No : HNSQ2020050063

LOT No : 001000
Date : 24 June 2025

Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
Humidification			
Humidity	50.0-90.0%	20.2 %	
Sterilization			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
Gas Injection			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
Exposure			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
Exhaust			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
Storage			
Storage Time	≥ 48 hour	≥ 48 hour	
Chemical Test			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
Biological Test			
EZ Test Biological Indicator	Green	Green	Pass


 (.....)
 Mr. Jutha Liam mukda
 Research and Development Engineer

			LOT	Color	BOX	TOTAL / PCS
		211098				
1	Z-LB-1.6-10.0-4	STERILE SURGICAL STEEL LABRET	LOT NO 01000			25
2	Z-LB-1.6-12.0-4	STERILE SURGICAL STEEL LABRET	LOT NO 01000			25
		211119				
1	Z-INMLBP289-1.2-6.0	STERILE SURGICAL STEEL INTERNALLY THREADED MICRO LABRET - HONEY BEE	LOT NO 01000			2
2	Z-MBLSCCM88-1.2-6.0-3-CR	STERILE SURGICAL STEEL HELIX MICRO BARBELL WITH JEWEL CHARM - STAR	LOT NO 01000			1
3	Z-MLBPE-1.2-6.0-2.5-WH	STERILE SURGICAL STEEL MICRO LABRET WITH SYNTHETIC PEARL	LOT NO 01000			1
		2542280				
1	Z-GPBNDJ-1.6-12.0-5/8-CR	STERILE GOLD PVD COATED SURGICAL STEEL DOUBLE PREMIUM CRYSTALS JEWEL Navel BANANA	LOT NO 01000			70
		2544418				
14	PT-Z-TINMLBJ-1.2-8.0-3-HP-CR	STERILE TITANIUM INTERNALLY THREADED JEWEL MICRO LABRET FOR PTIT ROC	LOT NO 01000			100
		2545211				
7	TDI-Z-TINMLB-1.2-8.0-3	STERILE TITANIUM INTERNALLY THREADED MICRO LABRET FOR TDI	LOT NO 01000			300
10	TDI-Z-TINMLBJ-1.2-8.0-3-CR	STERILE TITANIUM INTERNALLY THREADED JEWEL MICRO LABRET FOR TDI	LOT NO 01000			300
		2545984				
4	THD-Z-BKBL-1.6-22.0-4	STERILE BLACK PVD COATED SURGICAL STEEL BARBELL FOR THENDRO	LOT NO 01000			100
5	THD-Z-BKBL-1.6-34.0-4	STERILE BLACK PVD COATED SURGICAL STEEL BARBELL FOR THENDRO	LOT NO 01000			200
38	THD-Z-JNO-1.0-6.5-2.35-SA	STERILE SURGICAL STEEL JEWEL CURVED NOSE STUD FOR THENDRO	LOT NO 01000			100
		2546432				
2	Z-GPSCES776	STERILE GOLD PVD COATED SURGICAL STEEL EAR STUDS PAIR	LOT NO 01000			20
5	Z-PEBSTUDS-10-WH	STERILE PAIR OF SYNTHETIC PEARL BALL EAR STUDS	LOT NO 01000			10
		2546507				
10	PT-Z-TINMLB-1.2-8.0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO LABRET FOR PTIT ROC	LOT NO 01000			100
		2546749				
4	THD-Z-BKBN-1.6-8.0-4	STERILE BLACK PVD COATED SURGICAL STEEL CURVED BARBELL FOR THENDRO	LOT NO 01000			100
		2546971				
1	SV-STERILIZE	SERVICE STERILE	LOT NO 01000			20
						1474

Sterilization:
Temperature: 45.0℃
Humidity: 20.2RH
Pressures: -3.4KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 13:46

Add EO gas
Temperature: 44.8℃
Humidity: 20.2RH
Pressures: -4.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 13:46

Evacuation
Temperature: 44.4℃
Humidity: 20.1RH
Pressures: -0.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 13:41

Keeping Temp
Temperature: 43.5℃
Humidity: 20.1RH
Pressures: -1.6KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 12:41

Keep pressure qualified
Temperature: 42.1℃
Humidity: 20.1RH
Pressures: -58.7KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 12:39

Evacuation
Temperature: 41.7℃
Humidity: 20.1RH
Pressures: -3.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 12:28

Heating-up:
Temperature: 23.2℃
Humidity: 20.1RH
Pressures: -3.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 11:50

Setting:
Temperature: 45.0℃
Humidity: 55.0RH
Pressures: -60.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 11:50

2416125

END LOT 1000

Air replacement:
Temperature: 44.8℃
Humidity: 20.1RH
Pressures: -1.6KPa
Sterilization time: 0Hour
Air replace rounds: 1
round time: 0min(s)
2025-6-24 18:14

Air replacement:
Temperature: 45.9℃
Humidity: 20.1RH
Pressures: -1.6KPa
Sterilization time: 0Hour
Air replace rounds: 2
round time: 0min(s)
2025-6-24 17:54

Air replacement:
Temperature: 45.0℃
Humidity: 20.2RH
Pressures: -1.6KPa
Sterilization time: 0Hour
Air replace rounds: 3
round time: 0min(s)
2025-6-24 17:33

Air replacement:
Temperature: 45.6℃
Humidity: 20.1RH
Pressures: -1.6KPa
Sterilization time: 0Hour
Air replace rounds: 4
round time: 0min(s)
2025-6-24 17:13

Air replacement:
Temperature: 45.5℃
Humidity: 20.1RH
Pressures: -1.6KPa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 16:52

Sterilization:
Temperature: 44.4℃
Humidity: 20.2RH
Pressures: -0.0KPa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 16:47

Sterilization:
Temperature: 44.8℃
Humidity: 20.1RH
Pressures: -4.1KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-6-24 13:48

