

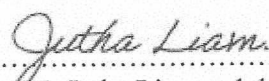
Report Test

M/C : Henan Sanqiang Medical Equipment Liability Co. Ltd
 Model : SQ-H220
 Serial No : HNSQ2020050063

LOT No : 001040
 Date : 25 August 2025

Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
Humidification			
Humidity	50.0-90.0%	20.2 %	
Sterilization			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
Gas Injection			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
Exposure			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
Exhaust			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
Storage			
Storage Time	≥ 48 hour	≥ 48 hour	
Chemical Test			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
Biological Test			
EZ Test Biological Indicator	Green	Green	Pass


 (.....)
 Mr. Jutha Liammukda
 Research and Development Engineer

			LOT	COLOR	BOX	TOTAL / PCS
		212901				
24	Z-LJNO-1.0-5.5-1.8-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELLED LARGE LEFT CURVE NOSE STUD	LOT NO 01040			80
		213279				
1	Z-TMLB-PIN-1.2-10.0-4B-HP	STERILE TITANIUM MICRO LABRET PIN	LOT NO 01040			50
2	Z-TMLB-PIN-1.2-6.0-4B-HP	STERILE TITANIUM MICRO LABRET PIN	LOT NO 01040			50
3	Z-TMLB-PIN-1.2-7.0-4B-HP	STERILE TITANIUM MICRO LABRET PIN	LOT NO 01040			50
4	Z-TMLB-PIN-1.2-8.0-4B-HP	STERILE TITANIUM MICRO LABRET PIN	LOT NO 01040			70
5	Z-TMLB-PIN-1.2-9.0-4B-HP	STERILE TITANIUM MICRO LABRET PIN	LOT NO 01040			50
		213296				
1	Z-GPMSL-1.2-8.0-3	STERILE GOLD PVD COATED SURGICAL STEEL MICRO BARBELL	LOT NO 01040			75
		213269				
8	Z-GPTBNH-1.6-10.0-4/8	STERILE GOLD PVD COATED TITANIUM NAVEL BANANA	LOT NO 01040			3
		213302				
1	Z-GPMMLB109-1.0-8.0-CR	STERILE GOLD PVD COATED SURGICAL STEEL INTERNALLY THREADED JEWELLED MICRO LABRET	LOT NO 01040			5
2	Z-GPMMLB123-1.2-8.0-CR	STERILE GOLD PVD COATED SURGICAL STEEL INTERNALLY THREADED JEWELLED MICRO LABRET - CRESCENT 3 GEMS	LOT NO 01040			5
3	Z-LJNO-1.0-5.0-2.35-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELLED LARGE LEFT CURVE NOSE STUD	LOT NO 01040			3
		213306				
1	Z-BNDJ-1.6-12.0-5/8-CR	STERILE SURGICAL STEEL DOUBLE JEWELLED NAVEL BANANA	LOT NO 01040			12
		213299				
3	Z-PCSH-1.0-8.0	STERILE SURGICAL STEEL HINGED SEGMENT RING	LOT NO 01040			3
4	Z-IMLB168-1.2-7.0-CR	STERILE SURGICAL STEEL INTERNALLY THREADED JEWELLED MICRO LABRET	LOT NO 01040			7
5	Z-IMLB23-1.0-6.0-CR	STERILE SURGICAL STEEL INTERNALLY THREADED JEWELLED MICRO LABRET	LOT NO 01040			5
6	Z-INO-0.6-6.5-2-CR	STERILE SURGICAL STEEL FLAT STONE PREMIUM CRYSTAL JEWELLED CURVED NOSE STUD	LOT NO 01040			3
7	Z-XBLU-8.1-6-10.0-5/8-BK-CR	STERILE BIOFLEX JEWELLED NAVEL BANANA WITH JEWELLED BLACK PVD BALL	LOT NO 01040			1
8	Z-XINLB11-1.2-8.0-3-CL-CR	STERILE BIOFLEX INTERNAL LABRET WITH JEWELLED SURGICAL STEEL ATTACHMENT	LOT NO 01040			15
9	Z-XINLB10-1.0-8.0-2-CL-CR	STERILE BIOFLEX INTERNAL MICRO LABRET WITH JEWELLED DISC	LOT NO 01040			5
10	Z-XINLBJS13-1.0-7.0-2-CL-CR	STERILE BIOFLEX INTERNAL LABRET WITH SILVER JEWELLED ATTACHMENT	LOT NO 01040			7
		2545894				
28	THD-Z-GPTNSCJ3-1.2-CR	STERILE GOLD PVD COATED TITANIUM MICRO JEWELLED ATTACHMENT FOR 1.2MM INTERNALLY THREADED PINS	LOT NO 01040			200
74	THD-Z-TMCCR-1.2-12.0-3	STERILE TITANIUM INTERNALLY THREADED MICRO CIRCULAR BARBELL	LOT NO 01040			100
		2553780				
11	PT-Z-TBL-1.6-20.0-5	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 01040			200
21	PT-Z-XBL-1.1-6-40.0-5-CL	STERILE BIOFLEX BARBELL WITH STEEL BALLS FOR PTIT ROC	LOT NO 01040			50
24	PT-Z-XBLB-1.1-2-12.0-3-4B-CL	STERILE BIOFLEX MICRO LABRET WITH STEEL BALL FOR PTIT ROC	LOT NO 01040			50
26	PT-Z10-TMRL-1.2-10.0-3	STERILE TITANIUM MICRO CURVED BARBELL FOR PTIT ROC	LOT NO 01040		SET/25	250
27	PT-Z10-TMLB-1.2-8.0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 01040		SET/125	1250
		2553940				
1	Z-GPTBNDJ-1.6-12.0-5/8-CR	STERILE GOLD PVD COATED TITANIUM DOUBLE PREMIUM CRYSTAL JEWELLED NAVEL BANANA	LOT NO 01040			15
2	Z-TBNJ-CH-1.6-10.0-5/8-HP-CR	STERILE TITANIUM DOUBLE VALUE JEWELLED NAVEL BANANA	LOT NO 01040			40
		2551570				
1	SV-STERILIZE	SERVICE STERILE	LOT NO 01040			1
		2553173				
1	SV-STERILIZE	SERVICE STERILE	LOT NO 01040			10
		2553438				
1	SV-STERILIZE	SERVICE STERILE	LOT NO 01040			15
		2553792				
2	SV-STERILIZE	SERVICE STERILE	LOT NO 01040			10
						2490

Sterilization:
Temperature: 45.0°C
Humidity: 20.2RH
Pressures: -0.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 13:41

Add EO gas
Temperature: 44.8°C
Humidity: 20.1RH
Pressures: -4.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 13:41

Evacuation
Temperature: 44.2°C
Humidity: 20.1RH
Pressures: -0.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 13:36

Keeping Temp
Temperature: 43.1°C
Humidity: 20.1RH
Pressures: -1.6KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 12:36

Keep pressure qualified
Temperature: 41.8°C
Humidity: 20.2RH
Pressures: -58.8KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 12:35

Evacuation
Temperature: 41.6°C
Humidity: 20.2RH
Pressures: -3.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 12:23

Heating-up:
Temperature: 23.8°C
Humidity: 20.1RH
Pressures: -3.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 11:46

Setting:
Temperature: 45.0°C
Humidity: 55.0RH
Pressures: -60.0KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 11:45

2518125
LOT 1040

END

Air replacement:
Temperature: 45.8°C
Humidity: 20.1RH
Pressures: -1.6KPa
Sterilization time: 0Hour
Air replace rounds: 1
round time: 0min(s)
2025-8-25 18:8

Air replacement:
Temperature: 45.3°C
Humidity: 20.1RH
Pressures: -1.6KPa
Sterilization time: 0Hour
Air replace rounds: 2
round time: 0min(s)
2025-8-25 17:48

Air replacement:
Temperature: 45.7°C
Humidity: 20.1RH
Pressures: -1.6KPa
Sterilization time: 0Hour
Air replace rounds: 3
round time: 0min(s)
2025-8-25 17:20

Air replacement:
Temperature: 46.0°C
Humidity: 20.2RH
Pressures: -1.6KPa
Sterilization time: 0Hour
Air replace rounds: 4
round time: 0min(s)
2025-8-25 17:8

Air replacement:
Temperature: 45.5°C
Humidity: 20.1RH
Pressures: -1.6KPa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 16:47

Sterilization:
Temperature: 45.2°C
Humidity: 20.1RH
Pressures: -0.5KPa
Sterilization time: 0Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 16:42

Sterilization:
Temperature: 44.7°C
Humidity: 20.1RH
Pressures: -4.1KPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2025-8-25 13:43