

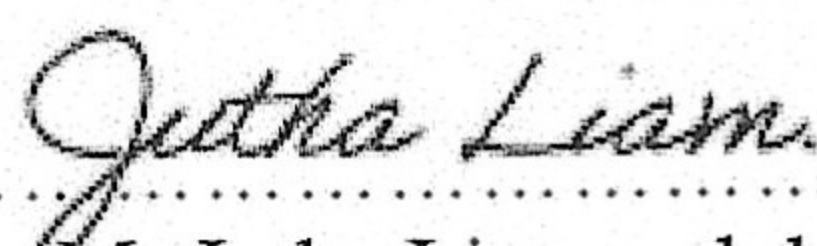
## Report Test

M/C : Henan Sanqiang Medical Equipment Liability Co. Ltd  
 Model : SQ-H220  
 Serial No : HNSQ2020050063

LOT No : 001058  
 Date : 19 September 2025

### Determination of Sterilization Operation Specification

| Item Detail                                                           | Specification | Setting   | Result |
|-----------------------------------------------------------------------|---------------|-----------|--------|
| Chamber Setting Temperature                                           | 40.0-55.5 °C  | 45.0 °C   |        |
| Vacuum pressure                                                       | 40-60 kPa     | 60 kPa    |        |
| Vacuum Time                                                           | 50-60 min     | 60 min    |        |
| <b>Humidification</b>                                                 |               |           |        |
| Humidity                                                              | 50.0-90.0%    | 20.2 %    |        |
| <b>Sterilization</b>                                                  |               |           |        |
| Sterilization Temperature                                             | 40.0-55.5 °C  | 43.5 °C   |        |
| <b>Gas Injection</b>                                                  |               |           |        |
| Final pressure                                                        | 0-10 kPa      | 4 kPa     |        |
| Final temperature of vaporizer                                        | 40.0-55.5 °C  | 44.8 °C   |        |
| <b>Exposure</b>                                                       |               |           |        |
| Temperature at start point of Exposure time period                    | 40.0-55.5 °C  | 43.3 °C   |        |
| Temperature at finish point of Exposure time period                   | 40.0-55.5 °C  | 46.0 °C   |        |
| Exposure Time                                                         | 120-240 min   | 180 min   |        |
| <b>Exhaust</b>                                                        |               |           |        |
| Exhaust Time                                                          | 60-80 min     | 75 min    |        |
| First Vacuum pressure                                                 | 0-2 kPa       | 1.5 kPa   |        |
| Second Vacuum pressure                                                | 0-2 kPa       | 1.5 kPa   |        |
| Third Vacuum pressure                                                 | 0-2 kPa       | 1.5 kPa   |        |
| Fourth Vacuum pressure                                                | 0-2 kPa       | 1.5 kPa   |        |
| Fifth Vacuum pressure                                                 | 0-2 kPa       | 1.5 kPa   |        |
| <b>Storage</b>                                                        |               |           |        |
| Storage Time                                                          | ≥ 48 hour     | ≥ 48 hour |        |
| <b>Chemical Test</b>                                                  |               |           |        |
| EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours | Blue          | Red       | Pass   |
| <b>Biological Test</b>                                                |               |           |        |
| EZ Test Biological Indicator                                          | Green         | Green     | Pass   |

  
 (.....)  
 Mr. Jutha Liam mukda  
 Research and Development Engineer

|    |                           |                                                                        | LOT         | Color | BOX    | TOTAL / PCS |
|----|---------------------------|------------------------------------------------------------------------|-------------|-------|--------|-------------|
|    |                           | 214126                                                                 |             |       |        |             |
| 2  | Z-TMLBJ-1.2-9.0-3-HP-CR   | STERILE TITANIUM PREMIUM CRYSTAL JEWELLED MICRO LABRET                 | LOT No 1058 |       |        | 20          |
|    |                           | 214235                                                                 |             |       |        |             |
| 51 | Z-MBN-1.2-10.0-3          | STERILE SURGICAL STEEL CURVED MICRO BARBELL                            | LOT No 1058 |       |        | 60          |
| 52 | Z-SES1-2.5-CR             | STERILE SURGICAL STEEL ROUND PRONG SET JEWELLED EAR STUDS PAIR         | LOT No 1058 |       |        | 50          |
| 53 | Z-SES1-3-CR               | STERILE SURGICAL STEEL ROUND PRONG SET JEWELLED EAR STUDS PAIR         | LOT No 1058 |       |        | 100         |
|    |                           | 214249                                                                 |             |       |        |             |
| 1  | Z-TBL-PIN-1.6-18.0-HP     | STERILE TITANIUM BARBELL PIN                                           | LOT No 1058 |       |        | 20          |
| 2  | Z-TBN-PIN-1.6-10.0-HP     | STERILE TITANIUM CURVED BARBELL PIN                                    | LOT No 1058 |       |        | 30          |
|    |                           | 214252                                                                 |             |       |        |             |
| 5  | ZDPU-1.5                  | STERILE DERMAL PUNCH                                                   | LOT No 1058 |       |        | 15          |
|    |                           | 214263                                                                 |             |       |        |             |
| 1  | Z-GPSCJNO12-0.8-8.5-CR    | STERILE GOLD PVD COATED SURGICAL STEEL CURVED PRONG SET JEWELLED TRINI | LOT No 1058 |       |        | 2           |
| 2  | Z-GPSCJNO53-0.8-8.5-CR    | STERILE GOLD PVD COATED SURGICAL STEEL CURVED JEWELLED NOSE STUD - FLO | LOT No 1058 |       |        | 2           |
| 3  | Z-GPSCJNO58-0.8-8.5-CR    | STERILE GOLD PVD COATED SURGICAL STEEL CURVED JEWELLED NOSE STUD - FLO | LOT No 1058 |       |        | 2           |
|    |                           | 214266                                                                 |             |       |        |             |
| 2  | Z-MBLDJ-1.2-8.0-3-CR      | STERILE SURGICAL STEEL DOUBLE JEWELLED MICRO BARBELL                   | LOT No 1058 |       |        | 100         |
|    |                           | 2556895                                                                |             |       |        |             |
| 5  | PT-Z-GPTMBL-1.2-10.0-3    | STERILE GOLD PVD COATED TITANIUM MICRO BARBELL                         | LOT No 1058 |       |        | 100         |
| 6  | PT-Z-GPTMBL-1.2-8.0-3     | STERILE GOLD PVD COATED TITANIUM MICRO BARBELL                         | LOT No 1058 |       |        | 100         |
| 8  | PT-Z-GPTMLB-1.2-8.0-3     | STERILE GOLD PVD COATED TITANIUM MICRO LABRET                          | LOT No 1058 |       |        | 100         |
| 16 | PT-Z-TJNO-0.8-8.5-2.35-CR | STERILE TITANIUM JEWELLED CURVED NOSE STUD FOR PTIT ROC                | LOT No 1058 |       |        | 250         |
| 17 | PT-Z-TMBN-1.2-8.0-3       | STERILE TITANIUM MICRO CURVED BARBELL FOR PTIT ROC                     | LOT No 1058 |       |        | 100         |
| 18 | PT-Z-TMCB-1.2-8.0-3       | STERILE TITANIUM MICRO CIRCULAR BARBELL FOR PTIT ROC                   | LOT No 1058 |       |        | 100         |
| 19 | PT-Z-TMLB-1.2-10.0-3      | STERILE TITANIUM MICRO LABRET FOR PTIT ROC                             | LOT No 1058 |       |        | 300         |
| 20 | PT-Z-TMLB-1.2-12.0-3      | STERILE TITANIUM MICRO LABRET FOR PTIT ROC                             | LOT No 1058 |       |        | 100         |
| 21 | PT-Z-TMLB-1.2-8.0-3       | STERILE TITANIUM MICRO LABRET FOR PTIT ROC                             | LOT No 1058 |       |        | 100         |
| 22 | PT-Z-TMLB-1.2-8.0-3       | STERILE TITANIUM MICRO LABRET FOR PTIT ROC                             | LOT No 1058 |       |        | 200         |
| 25 | PT-Z10-TBL-1.6-22.0-5     | STERILIZED TITANIUM BARBELL FOR PTIT ROC                               | LOT No 1058 |       | SET/25 | 250         |
| 28 | PT-Z10-TMLBJ-1.2-8.0-3-CR | STERILE TITANIUM MICRO LABRET WITH JEWELLED BALL FOR PTIT ROC          | LOT No 1058 |       | SET/25 | 250         |
|    |                           | 2556831                                                                |             |       |        |             |
| 1  | TDI-Z-TMLB-1.2-10.0-3     | STERILE TITANIUM MICRO LABRET FOR TDI                                  | LOT No 1058 |       |        | 300         |
|    |                           |                                                                        |             |       |        | 2651        |

Sterilization:  
Temperature: 43.9t  
Humidity: 20.1RH  
Pressures: -3.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 14:39

Add ED gas  
Temperature: 43.9t  
Humidity: 20.2RH  
Pressures: -4.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 14:39

Evacuation  
Temperature: 44.2t  
Humidity: 20.1RH  
Pressures: -0.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 13:46

Keeping Temp  
Temperature: 42.7t  
Humidity: 20.1RH  
Pressures: -1.6KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 12:46

Keep pressure qualified  
Temperature: 41.4t  
Humidity: 20.1RH  
Pressures: -56.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 12:44

Evacuation  
Temperature: 41.5t  
Humidity: 20.2RH  
Pressures: -3.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 12:32

Heating-up:  
Temperature: 22.2t  
Humidity: 20.1RH  
Pressures: -3.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 11:50

Setting:  
Temperature: 45.0t  
Humidity: 55.0RH  
Pressures: -60.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 11:50

19/9/25  
ENC Lot 1058

Air replacement:  
Temperature: 44.6t  
Humidity: 20.1RH  
Pressures: -1.6KPa  
Sterilization time: 0Hour  
Air replace rounds: 1  
round time: 0min(s)  
2025-9-19 19:7

Air replacement:  
Temperature: 44.8t  
Humidity: 20.1RH  
Pressures: -1.6KPa  
Sterilization time: 0Hour  
Air replace rounds: 2  
round time: 0min(s)  
2025-9-19 13:46

Air replacement:  
Temperature: 45.1t  
Humidity: 20.1RH  
Pressures: -1.6KPa  
Sterilization time: 0Hour  
Air replace rounds: 3  
round time: 0min(s)  
2025-9-19 13:26

Air replacement:  
Temperature: 46.3t  
Humidity: 20.1RH  
Pressures: -1.6KPa  
Sterilization time: 0Hour  
Air replace rounds: 4  
round time: 0min(s)  
2025-9-19 16:6

Air replacement:  
Temperature: 45.7t  
Humidity: 20.1RH  
Pressures: -1.6KPa  
Sterilization time: 0Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 17:46

Sterilization:  
Temperature: 45.0t  
Humidity: 20.2RH  
Pressures: -2.8KPa  
Sterilization time: 0Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 17:40

Sterilization:  
Temperature: 44.0t  
Humidity: 20.1RH  
Pressures: -4.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2025-9-19 14:41