

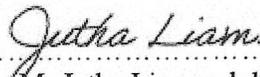
## Report Test

M/C : Henan Sanqiang Medical Equipment Liability Co. Ltd  
 Model : SQ-H220  
 Serial No : HNSQ2020050063

LOT No : 001056  
 Date : 17 September 2025

### Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
<b>Humidification</b>			
Humidity	50.0-90.0%	20.2 %	
<b>Sterilization</b>			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
<b>Gas Injection</b>			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
<b>Exposure</b>			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
<b>Exhaust</b>			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
<b>Storage</b>			
Storage Time	≥ 48 hour	≥ 48 hour	
<b>Chemical Test</b>			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
<b>Biological Test</b>			
EZ Test Biological Indicator	Green	Green	Pass

  
 (.....)  
 Mr. Jutha Liammukda  
 Research and Development Engineer



FM-ISO0-009 Rev.02 (05/11/2020)

			LOT	Color	BOX	TOTAL / PCS
		214109				
147	Z-SOBYC10-1.6-CR	STERILE SURGICAL STEEL JEWELED BELLY CLICKER	LOT NO 1056			10
148	Z-SCSH225-1.2-10.0-CR-AQ	STERILE SURGICAL STEEL JEWELED MULTI PURPOSE CLICKER	LOT NO 1056			10
		214124				
8	Z-JMB-1.2-2.5-AM	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELED MICRO BALL	LOT NO 1056			6
9	Z-JMB-1.2-2.5-AQ	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELED MICRO BALL	LOT NO 1056			9
10	Z-JMB-1.2-2.5-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELED MICRO BALL	LOT NO 1056			10
11	Z-JMB-1.2-3-AM	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELED MICRO BALL	LOT NO 1056			5
12	Z-JMB-1.2-3-JE	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELED MICRO BALL	LOT NO 1056			6
		214131				
1	Z-BL-1.6-22.0-5	STERILE SURGICAL STEEL BARBELL	LOT NO 1056			100
2	Z-MBL-1.2-8.0-3	STERILE SURGICAL STEEL MICRO BARBELL	LOT NO 1056			200
		214198				
1	Z-MBPACK6-1.2-3	PACK OF 6 STERILE SURGICAL STEEL MICRO BALLS	LOT NO 1056			2
2	Z-MLB-1.2-10.0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 1056			10
3	Z-MLB-1.2-8.0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 1056			10
		2551120				
27	THD-Z-SCJNO12-1.0-6.5-CR	STERILE SURGICAL STEEL CURVED PRONG SET JEWELED TRINITY NOSE STUD	LOT NO 1056			150
		2556179				
8	PT-Z-TINMBN-1.2-10.0-3-HP	STERILE TITANIUM INTERNALLY THREADED CURVED MICRO BARBELL	LOT NO 1056			50
18	PT-Z10-TMCB-1.2-10.0-3	STERILE TITANIUM MICRO CIRCULAR BARBELL FOR PITT ROC	LOT NO 1056		SET/50	500
		2566644				
1	TDI-Z-TBCSH-1.2-8.0	STERILE TITANIUM HINGED SEGMENT RING FOR TDI	LOT NO 1056			300
						1378

Sterilization:  
Temperature: 45.7t  
Humidity:20.1RH  
Pressures:-4.1KPa  
Sterilization time:3Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 13:54

Add EO gas  
Temperature: 45.6t  
Humidity:20.1RH  
Pressures:-4.0KPa  
Sterilization time:3Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 13:54

Evacuation  
Temperature: 44.2t  
Humidity:20.1RH  
Pressures:-0.0KPa  
Sterilization time:3Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 13:44

Keeping Temp  
Temperature: 43.0t  
Humidity:20.1RH  
Pressures:-1.6KPa  
Sterilization time:3Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 12:44

Keep pressure qualified  
Temperature: 41.5t  
Humidity:20.1RH  
Pressures:-58.7KPa  
Sterilization time:3Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 12:43

Evacuation  
Temperature: 41.5t  
Humidity:20.1RH  
Pressures:-3.0KPa  
Sterilization time:3Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 12:31

Heating-up:  
Temperature: 22.2t  
Humidity:20.2RH  
Pressures:-3.0KPa  
Sterilization time:3Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 11:50

Setting:  
Temperature: 45.0t  
Humidity:55.0RH  
Pressures:-60.0KPa  
Sterilization time:3Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 11:50

19/09/25

END  
Lot 1005b

Air replacement:  
Temperature: 45.3t  
Humidity:20.1RH  
Pressures:-1.6KPa  
Sterilization time:0Hour  
Air replace rounds:1  
round time:0min(s)  
2025-9-17 10:21

Air replacement:  
Temperature: 45.0t  
Humidity:20.1RH  
Pressures:-.5KPa  
Sterilization time:0Hour  
Air replace rounds:2  
round time:0min(s)  
2025-9-17 10:1

Air replacement:  
Temperature: 45.6t  
Humidity:20.1RH  
Pressures:-1.6KPa  
Sterilization time:0Hour  
Air replace rounds:3  
round time:0min(s)  
2025-9-17 11:41

Air replacement:  
Temperature: 45.6t  
Humidity:20.1RH  
Pressures:-1.6KPa  
Sterilization time:0Hour  
Air replace rounds:4  
round time:0min(s)  
2025-9-17 17:21

Air replacement:  
Temperature: 45.0t  
Humidity:20.1RH  
Pressures:-1.6KPa  
Sterilization time:0Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 17:0

Sterilization:  
Temperature: 45.0t  
Humidity:20.1RH  
Pressures:-2.3KPa  
Sterilization time:0Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 16:55

Sterilization:  
Temperature: 45.5t  
Humidity:20.1RH  
Pressures:-5.2KPa  
Sterilization time:3Hour  
Air replace rounds:5  
round time:15min(s)  
2025-9-17 13:56

ISO 9001



ISO 14001



ISO 45001



SEDEX



CSR-DIW



FM-ISO-009 Rev.02 (05/11/2020)