

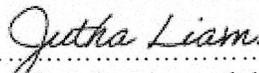
Report Test

M/C: Henan Sanqiang Medical Equipment Liability Co. Ltd
 Model: SQ-H220
 Serial No: HNSQ2020050063

LOT No: 001111
 Date: 11 December 2025

Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
Humidification			
Humidity	50.0-90.0%	20.2 %	
Sterilization			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
Gas Injection			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
Exposure			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
Exhaust			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
Storage			
Storage Time	≥ 48 hour	≥ 48 hour	
Chemical Test			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
Biological Test			
EZ Test Biological Indicator	Green	Green	Pass


 (.....)
 Mr. Jutha Liamukda
 Research and Development Engineer

			LOT	Color	BOX	TOTAL / PCS
			216622			
1	Z-GPBNJ-1.6-10.0-5-8-CR	STERILE GOLD PVD COATED SURGICAL STEEL PREMIUM CRYSTAL JEWELLED NAVEL BANANA	LOT NO 1111			50
			216624			
1	Z-JNO-1.0-6.5-2.35-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELLED CURVED NOSE STUD	LOT NO 1111			200
2	Z-TBNDJ-1.6-10.0-5/8-HP-CR	STERILE TITANIUM DOUBLE JEWELLED NAVEL BANANA	LOT NO 1111			50
			216706			
7	Z-TESI-4-JE	STERILE PAIR OF TITANIUM ROUND PRONG SET JEWELLED EAR STUDS	LOT NO 1111			5
			216707			
1	Z-TESI-5-JE	STERILE PAIR OF TITANIUM ROUND PRONG SET JEWELLED EAR STUDS	LOT NO 1111			5
			216725			
2	Z-BL-1.6-10.0-8	STERILE SURGICAL STEEL BARBELL	LOT NO 1111			1
3	Z-BL-1.6-12.0-3	STERILE SURGICAL STEEL BARBELL	LOT NO 1111			1
			216728			
1	Z-BL-1.6-18.0-5	STERILE SURGICAL STEEL BARBELL	LOT NO 1111			10
2	Z-BL-1.6-20.0-5	STERILE SURGICAL STEEL BARBELL	LOT NO 1111			10
3	Z-CNOJ90-0.6-6.5-1.7-CR	STERILE SURGICAL STEEL 90 DEGREE FLAT PREMIUM CRYSTAL JEWELLED NOSE STUD	LOT NO 1111			10
4	Z-MLB-1.2-8.0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 1111			20
			216729			
1	Z-ZGJNO90-1.0-6.5-2.35-CR	STERILE ZIRCON GOLD PVD COATED SURGICAL STEEL PREMIUM CRYSTAL JEWELLED 90 DEGREE NOSE STUD	LOT NO 1111			5
			216730			
4	Z-LJNO-1.0-6.5-1.8-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELLED LARGE LEFT CURVE NOSE STUD	LOT NO 1111			6
5	Z-ZGMLBJ-1.2-10.0-3-CR	STERILE ZIRCON GOLD PVD SURGICAL STEEL JEWELLED MICRO LABRET	LOT NO 1111			6
6	Z-ZGMLB-1.2-8.0-3-CR	STERILE ZIRCON GOLD PVD SURGICAL STEEL JEWELLED MICRO LABRET	LOT NO 1111			6
			216742			
21	Z-BNDJ-1.6-10.0-4/8-BZ	STERILE SURGICAL STEEL DOUBLE PREMIUM CRYSTALS JEWELLED NAVEL BANANA	LOT NO 1111			20
			216748			
1	Z-BC-1.6-8.0-4	STERILE SURGICAL STEEL BCR	LOT NO 1111			4
2	Z-GPSCS430-3	STERILE GOLD PVD COATED SURGICAL STEEL EAR STUDS PAIR - BALL 3MM	LOT NO 1111			1
3	Z-LB-1.6-6.0-3	STERILE SURGICAL STEEL LABRET	LOT NO 1111			4
4	Z-LB-1.0-7.0-3	STERILE SURGICAL STEEL LABRET	LOT NO 1111			4
5	Z-SBNDJ-CN-1.6-8.0-4/8-CR	STERILE SURGICAL STEEL DOUBLE VALUE CRYSTAL MINI NAVEL BANANA	LOT NO 1111			2
6	Z-SCS430-4	STERILE SURGICAL STEEL EAR STUDS PAIR - BALL 3MM	LOT NO 1111			3
7	Z-TBCX-3.0-10.0-5-HP	STERILE TITANIUM BALL CLOSURE RING WITH POP OUT BALL	LOT NO 1111			1
8	Z-TINMBL-1.2-8.0-3-HP	STERILE TITANIUM INTERNALLY THREADED CURVED MICRO BARBELL	LOT NO 1111			2
9	Z-TINICB-1.2-10.0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO CIRCULAR BARBELL	LOT NO 1111			2
10	Z-TIMBL-1.2-38.0-3-HP	STERILE TITANIUM MICRO BARBELL	LOT NO 1111			2
			216747			
1	Z-LB-1.6-8.0-4	STERILE SURGICAL STEEL LABRET	LOT NO 1111			4
2	Z-TINMBN-1.2-8.0-3-HP	STERILE TITANIUM INTERNALLY THREADED CURVED MICRO BARBELL	LOT NO 1111			2
			216749			
2	Z-INMLBJ-1.2-8.0-2-AM	STERILE SURGICAL STEEL INTERNALLY THREADED PREMIUM CRYSTAL JEWELLED MICRO LABRET	LOT NO 1111			1
3	Z-INMLSC45-1.2-8.0-3	STERILE SURGICAL STEEL TRAGUS MICRO BARBELL - SALAMANDER	LOT NO 1111			2
4	Z-INMLSCCM88-1.2-6.0-3-CR	STERILE SURGICAL STEEL HELIX MICRO BARBELL WITH JEWELLED CHARM - STAR	LOT NO 1111			2
5	Z-INMSC307-1.2	STERILE SURGICAL STEEL MICRO THREADED ATTACHMENT	LOT NO 1111			1
6	Z-TINMLBJD-1.0-8.0-2-HP-CR	STERILE TITANIUM INTERNALLY THREADED JEWELLED DISC MICRO LABRET	LOT NO 1111			2
7	Z-XINLBSJD-1.2-8.0-2.35-CL-CR	STERILE BIOFLEX INTERNAL LABRET WITH SURGICAL STEEL JEWELLED DISC	LOT NO 1111			2
			216780			
1	Z-TBL-1.6-16.0-5-HP	STERILE TITANIUM BARBELL	LOT NO 1111			76
			2563413			
3	THD-Z-BKINMLBJ-1.2-10.0-3-CR	STERILE BLACK PVD COATED SURGICAL STEEL INTERNALLY THREADED PREMIUM CRYSTAL JEWELLED MICRO LABRET	LOT NO 1111			100
19	THD-Z-GPMBLSJ-1.2-8.0-3-CR	STERILE GOLD PVD COATED SURGICAL STEEL ROUND PRONG SET JEWELLED TRAGUS MICRO BARBELL	LOT NO 1111			130
26	THD-Z-GPTINMB-1.2-3	STERILE GOLD PVD COATED TITANIUM MICRO BALL FOR 1.2MM INTERNALLY THREADED PINS	LOT NO 1111			200
			2563421			
3	THD-Z-14NOJ19-C-0.8-6.5-CR	STERILE 14K GOLD CURVED PRONG SET FLOWER NOSE STUD	LOT NO 1111			30
			2564378			
11	THD-Z-TJNO-PL-1.0-6.5-2.35-WOP	STERILE TITANIUM SYNTHETIC OPAL CURVED NOSE STUD	LOT NO 1111			100
			2565367			
24	PT-Z-XJNO2-0.8-6.0-2.35-CL-CR	STERILE BIOFLEX JEWELLED CURVED NOSE STUD FOR PTTIT ROC	LOT NO 1111			50
			2565712			
1	TDI-Z-TINDAID-1.6-1.5-4-CR	STERILE TITANIUM INTERNALLY THREADED DERMAL ANCHOR WITH JEWELLED DISC FOR TDI	LOT NO 1111			100
			2565853			
1	PT-Z-TMLBJ-1.2-8.0-3-CR	STERILE TITANIUM MICRO LABRET WITH JEWELLED BALL FOR PTTIT ROC	LOT NO 1111			100
			2562598			
1	SV-STERILIZE	SERVICE STERILE	LOT NO 1111			90
						1422

Sterilization:
Temperature: 45.5t
Humidity: 20.1RH
Pressures: -0.6kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 13:31

Add E0 gas
Temperature: 45.4t
Humidity: 20.1RH
Pressures: -4.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 15:01

Evacuation
Temperature: 44.3t
Humidity: 20.1RH
Pressures: -0.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 14:55

Keeping Temp
Temperature: 43.6t
Humidity: 20.1RH
Pressures: -1.6kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 13:55

Keep pressure qualified
Temperature: 42.1t
Humidity: 20.1RH
Pressures: -0.8kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 13:53

Evacuation
Temperature: 41.7t
Humidity: 20.1RH
Pressures: -0.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 13:41

Heating-up:
Temperature: 38.2t
Humidity: 20.1RH
Pressures: -0.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 13:15

Setting:
Temperature: 45.0t
Humidity: 55.0RH
Pressures: -50.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 13:15

END

11/12/25
Lot 4447

Air replacement:
Temperature: 45.2t
Humidity: 20.1RH
Pressures: -1.6kPa
Sterilization time: 0hour
Air replace rounds: 1
round time: 0min(s)
2025-12-11 13:58

Air replacement:
Temperature: 45.0t
Humidity: 20.1RH
Pressures: -1.6kPa
Sterilization time: 0hour
Air replace rounds: 2
round time: 0min(s)
2025-12-11 13:58

Air replacement:
Temperature: 45.3t
Humidity: 20.1RH
Pressures: -1.6kPa
Sterilization time: 0hour
Air replace rounds: 3
round time: 0min(s)
2025-12-11 13:18

Air replacement:
Temperature: 45.2t
Humidity: 20.1RH
Pressures: -1.6kPa
Sterilization time: 0hour
Air replace rounds: 4
round time: 0min(s)
2025-12-11 13:57

Air replacement:
Temperature: 45.5t
Humidity: 20.1RH
Pressures: -1.6kPa
Sterilization time: 0hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 13:07

Sterilization:
Temperature: 45.1t
Humidity: 20.1RH
Pressures: -0.6kPa
Sterilization time: 0hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 13:22

Sterilization:
Temperature: 45.5t
Humidity: 20.1RH
Pressures: -4.1kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2025-12-11 15:33