

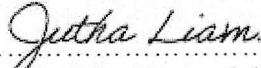
## Report Test

M/C: Henan Sanqiang Medical Equipment Liability Co. Ltd  
 Model: SQ-H220  
 Serial No: HNSQ2020050063

LOT No: 001123  
 Date: 07 January 2026

### Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
<b>Humidification</b>			
Humidity	50.0-90.0%	20.2 %	
<b>Sterilization</b>			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
<b>Gas Injection</b>			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
<b>Exposure</b>			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
<b>Exhaust</b>			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
<b>Storage</b>			
Storage Time	≥ 48 hour	≥ 48 hour	
<b>Chemical Test</b>			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
<b>Biological Test</b>			
EZ Test Biological Indicator	Green	Green	Pass

  
 (.....)  
 Mr. Jutha Liammukda  
 Research and Development Engineer

		LOT	Color	BOX	TOTAL / PCS
		217234			
1	Z-BNDJ-1.6-12.0-5/8-CR	STERILE SURGICAL STEEL DOUBLE JEWELLED NAVEL BANANA	LOT NO 001123		50
2	Z-SES1-4-CR	STERILE SURGICAL STEEL ROUND PRONG SET JEWELLED EAR STUDS PAIR	LOT NO 001123		100
4	Z-TMLBJ-1.2-8.0-3-HP-CR	STERILE TITANIUM MICRO LABRET WITH JEWELLED BALL	LOT NO 001123		50
6	Z-XBNJ-1.1.6-12.0-5/8-CL-CR	BIOFLEX JEWELLED NAVEL BANANA WITH SURGICAL STEEL BALL	LOT NO 001123		100
		217297			
1	Z-BL-1.6-12.7-4	STERILE SURGICAL STEEL BARBELL	LOT NO 001123		10
2	Z-BL-1.6-14.0-4	STERILE SURGICAL STEEL BARBELL	LOT NO 001123		5
3	Z-BL-1.6-16.0-4	STERILE SURGICAL STEEL BARBELL	LOT NO 001123		5
4	Z-BL-1.6-44.0-4	STERILE SURGICAL STEEL BARBELL	LOT NO 001123		5
5	Z-BNDJ-1.6-12.0-4/8-CR	STERILE SURGICAL STEEL DOUBLE PREMIUM CRYSTALS JEWELLED NAVEL BANANA	LOT NO 001123		10
6	Z-EX-2.5	STERILE SURGICAL STEEL EXPANDER	LOT NO 001123		2
7	Z-EX-4.0	STERILE SURGICAL STEEL EXPANDER	LOT NO 001123		2
8	Z-EX-6.0	STERILE SURGICAL STEEL EXPANDER	LOT NO 001123		2
9	Z-EX-9.0	STERILE SURGICAL STEEL EXPANDER	LOT NO 001123		2
10	Z-GPBNJ-1.6-12.0-5/8-CR	STERILE GOLD PVD COATED SURGICAL STEEL PREMIUM CRYSTAL JEWELLED NAVEL BANANA	LOT NO 001123		5
11	Z-GPJNO-0.8-6.5-1.8-CR	STERILE GOLD PVD COATED SURGICAL STEEL PREMIUM CRYSTAL JEWELLED CURVED NOSE STUD	LOT NO 001123		10
12	Z-JNO-0.5-6.5-1.8-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELLED CURVED NOSE STUD	LOT NO 001123		10
13	Z-LB-1.6-8.0-3	STERILE SURGICAL STEEL LABRET	LOT NO 001123		10
14	Z-MBN-1.2-10.0-3	STERILE SURGICAL STEEL CURVED MICRO BARBELL	LOT NO 001123		10
15	Z-MLB-1.0-8.0-2.5	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 001123		10
16	Z-MLB-1.2-8.0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 001123		10
17	Z-MLBJ-1.0-6.0-3-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELLED MICRO LABRET	LOT NO 001123		10
18	Z-MLBJ-1.0-8.0-3-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELLED MICRO LABRET	LOT NO 001123		10
19	Z-ZGBL-1.6-14.0-5	STERILE ZIRCON GOLD PVD COATED SURGICAL STEEL BARBELL	LOT NO 001123		2
20	Z-ZGBL-1.6-16.0-4	STERILE ZIRCON GOLD PVD COATED SURGICAL STEEL BARBELL	LOT NO 001123		2
21	Z-ZGMLB-1.2-8.0-3	STERILE ZIRCON GOLD PVD COATED SURGICAL STEEL MICRO LABRET	LOT NO 001123		10
22	Z-ZGMLB-1.2-8.0-3	STERILE ZIRCON GOLD PVD COATED SURGICAL STEEL MICRO LABRET	LOT NO 001123		10
23	Z-ZGMLBJ-1.2-6.0-3-CR	STERILE ZIRCON GOLD PVD COATED SURGICAL STEEL JEWELLED MICRO LABRET	LOT NO 001123		10
24	Z-ZGMLBJ-1.2-8.0-3-CR	STERILE ZIRCON GOLD PVD COATED SURGICAL STEEL JEWELLED MICRO LABRET	LOT NO 001123		10
		217327			
2	Z-TBN-PIN-1.6-10.0-HP	STERILE TITANIUM CURVED BARBELL PIN	LOT NO 001123		50
3	Z-TMLB-PIN-1.2-7.0-4B-HP	STERILE TITANIUM MICRO LABRET PIN	LOT NO 001123		100
4	Z-TMLB-PIN-1.2-8.0-4B-HP	STERILE TITANIUM MICRO LABRET PIN	LOT NO 001123		100
5	Z-TMLB-PIN-1.2-9.0-4B-HP	STERILE TITANIUM MICRO LABRET PIN	LOT NO 001123		100
		2562833			
2	PT-Z-GPTBNJ-1.6-12.0-5/8-CR	STERILE GOLD PVD COATED TITANIUM PREMIUM CRYSTAL JEWELLED NAVEL BANANA	LOT NO 001123		100
		2563413			
17	THD-Z-GPINSJ3L-1.6-5-WOP	STERILE GOLD PVD COATED SURGICAL STEEL SYNTHETIC OPAL ATTACHMENT FOR 1.6MM INTERNAL	LOT NO 001123		100
					1022



FM-ISO0-009 Rev.02 (05/11/2020)

Sterilization:  
Temperature: 45.5°C  
Humidity: 20.28%  
Pressures: -2.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 13:12

Add EO gas  
Temperature: 45.2°C  
Humidity: 20.18%  
Pressures: -4.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 13:11

Evacuation  
Temperature: 44.1°C  
Humidity: 20.18%  
Pressures: -0.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 13:4

Keeping Temp  
Temperature: 42.9°C  
Humidity: 20.18%  
Pressures: -1.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 12:4

Keep pressure qualified  
Temperature: 41.5°C  
Humidity: 20.18%  
Pressures: -0.3KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 12:2

Evacuation  
Temperature: 41.0°C  
Humidity: 20.18%  
Pressures: -0.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 11:50

Heating up  
Temperature: 21.5°C  
Humidity: 20.18%  
Pressures: -0.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 11:6

Setting  
Temperature: 45.0°C  
Humidity: 20.18%  
Pressures: -0.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 11:6

IND 7/01/26  
Lot 1123

Air replacement:  
Temperature: 45.8°C  
Humidity: 20.18%  
Pressures: -1.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 1  
round time: 15min(s)  
2026-1-7 11:39

Air replacement:  
Temperature: 45.2°C  
Humidity: 20.18%  
Pressures: -1.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 2  
round time: 15min(s)  
2026-1-7 11:19

Air replacement:  
Temperature: 46.1°C  
Humidity: 20.18%  
Pressures: -1.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 3  
round time: 15min(s)  
2026-1-7 10:59

Air replacement:  
Temperature: 45.9°C  
Humidity: 20.28%  
Pressures: -1.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 4  
round time: 15min(s)  
2026-1-7 10:39

Air replacement:  
Temperature: 45.3°C  
Humidity: 20.18%  
Pressures: -1.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 10:19

Sterilization:  
Temperature: 45.2°C  
Humidity: 20.28%  
Pressures: -0.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 10:12

Sterilization:  
Temperature: 45.4°C  
Humidity: 20.18%  
Pressures: -0.0KPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-7 11:14

