

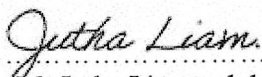
Report Test

M/C: Henan Sanqiang Medical Equipment Liability Co. Ltd
 Model: SQ-H220
 Serial No: HNSQ2020050063

LOT No: 001128
 Date: 14 January 2026

Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
Humidification			
Humidity	50.0-90.0%	20.2 %	
Sterilization			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
Gas Injection			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
Exposure			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
Exhaust			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
Storage			
Storage Time	≥ 48 hour	≥ 48 hour	
Chemical Test			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
Biological Test			
EZ Test Biological Indicator	Green	Green	Pass


 (.....)
 Mr. Jutha Liammukda
 Research and Development Engineer



FM-ISO0-009 Rev.02 (05/11/2020)

			LOT	Color	BOX	TOTAL / PCS
		217234				
6	Z-TPFINJD-1 0-8 0-3-HP-CR	STERILE TITANIUM JEWELED THREADLESS DISC LABRET	LOT NO 001128			40
		217589				
2	Z-TINDAJD-1 6 2.5-4-HP-CR	STERILE TITANIUM INTERNALLY THREADED DERMAL ANCHOR WITH JEWELED DISC	LOT NO 001128			25
3	Z-TMBN-1 2-10 0-3-HP	STERILE TITANIUM MICRO CURVED BARBELL	LOT NO 001128			75
		217593				
1	Z-SES1-2 5-CR	STERILE SURGICAL STEEL ROUND PRONG SET JEWELED EAR STUDS PAIR	LOT NO 001128			15
		217598				
2	Z-TINMLBJ-1 2-8 0-2-HP-CR	STERILE TITANIUM INTERNALLY THREADED JEWELED MICRO LABRET	LOT NO 001128			10
		217527				
21	Z-TB-1 6-3-HP	STERILE TITANIUM BALL	LOT NO 001128			8
22	Z-TB-1 6-4-HP	STERILE TITANIUM BALL	LOT NO 001128			5
23	Z-TB-1 6-5-HP	STERILE TITANIUM BALL	LOT NO 001128			6
24	Z-TBNDJ-1 6-10 0-5/8-HP-AQ	STERILE TITANIUM DOUBLE JEWELED NAVEL BANANA	LOT NO 001128			4
25	Z-TBNDJ-1 6-10 0-5/8-HP-SA	STERILE TITANIUM DOUBLE JEWELED NAVEL BANANA	LOT NO 001128			2
26	Z-TINMLB-PIN-1 2-6 0-4B-HP	STERILE TITANIUM INTERNALLY THREADED MICRO LABRET PIN	LOT NO 001128			8
		2565367				
10	PT-Z-GPTMLBJ-1 2-8 0-3-CR	STERILE GOLD PVD COATED TITANIUM PREMIUM CRYSTAL JEWELED MICRO LABRET	LOT NO 001128			300
		2566800				
7	PT-Z-GPTMLBJ-1 2-8 0-3-CR	STERILE GOLD PVD COATED TITANIUM PREMIUM CRYSTAL JEWELED MICRO LABRET	LOT NO 001128			200
		2567330				
1	Z-TPFINUB-1 0-8 0-3-HP-CR	STERILE TITANIUM JEWELED THREADLESS LABRET	LOT NO 001128			20
		2568458				
39	PT-Z10-TMBN-1 2-10 0-3	STERILE TITANIUM MICRO CURVED BARBELL FOR PTIT ROC	LOT NO 001128		SET50	500
43	PT-Z10-TMLB-1 2-9 0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 001128		SET100	1000
46	PT-Z10-TMLBJ-1 2-9 0-2.5-CR	STERILE TITANIUM MICRO LABRET WITH JEWELED BALL FOR PTIT ROC	LOT NO 001128		SET50	500
		2568899				
1	Z-TINMLBJ-1 2-10 0-2-HP-CR	STERILE TITANIUM INTERNALLY THREADED JEWELED MICRO LABRET	LOT NO 001128			10
		2567982				
2	SV-STERILIZE	SERVICE STERILE	LOT NO 001128			21
		2564106				
1	ZNEEDLES-1 0-48 0	SLBJ STOCK ONLY	LOT NO 001128			100
2	ZNEEDLES-1 2-48 0	SLBJ STOCK ONLY	LOT NO 001128			100
3	ZNEEDLES-1 6-48 0	SLBJ STOCK ONLY	LOT NO 001128			100
		2564111				
1	ZNEEDLES-1 0-48 0	SLBJ STOCK ONLY	LOT NO 001128			100
2	ZNEEDLES-1 2-48 0	SLBJ STOCK ONLY	LOT NO 001128			100
3	ZNEEDLES-1 6-48 0	SLBJ STOCK ONLY	LOT NO 001128			100
		2564418				
1	ZNEEDLES-1 2-48 0	SLBJ STOCK ONLY	LOT NO 001128			300
		2567728				
1	ZNEEDLES-1 0-48 0	SLBJ STOCK ONLY	LOT NO 001128			200
2	ZNEEDLES-1 2-48 0	SLBJ STOCK ONLY	LOT NO 001128			200
3	ZNEEDLES-1 6-48 0	SLBJ STOCK ONLY	LOT NO 001128			200
		2565907				
2	ZNEEDLES-1 2-48 0	SLBJ STOCK ONLY	LOT NO 001128			100
		2566680				
1	ZNEEDLES-1 0-48 0	SLBJ STOCK ONLY	LOT NO 001128			100
2	ZNEEDLES-1 2-48 0	SLBJ STOCK ONLY	LOT NO 001128			100
3	ZNEEDLES-1 6-48 0	SLBJ STOCK ONLY	LOT NO 001128			100
		2567418				
1	ZNEEDLES-1 0-48 0	SLBJ STOCK ONLY	LOT NO 001128			500
2	ZNEEDLES-1 2-48 0	SLBJ STOCK ONLY	LOT NO 001128			400
3	ZNEEDLES-1 6-48 0	SLBJ STOCK ONLY	LOT NO 001128			200
						5749

Sterilization:
Temperature: 43.7°C
Humidity: 20.1RH
Pressure: -2.16Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-14 13:46

Add (G) gas
Temperature: 43.5°C
Humidity: 20.1RH
Pressure: -4.09Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-14 13:48

Evacuation
Temperature: 44.1°C
Humidity: 20.1RH
Pressure: -0.00Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-14 13:42

Keeping Temp
Temperature: 42.7°C
Humidity: 20.1RH
Pressure: -1.00Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-14 12:42

Keep pressure qualified
Temperature: 41.4°C
Humidity: 20.1RH
Pressure: -59.18Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-14 12:40

Evacuation
Temperature: 41.5°C
Humidity: 20.1RH
Pressure: -0.00Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-14 12:38

Heating-up
Temperature: 21.5°C
Humidity: 20.1RH
Pressure: -3.00Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-14 11:45

Setting
Temperature: 45.0°C
Humidity: 55.0RH
Pressure: -60.00Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-14 11:45

END

14:09 Lab
Lab 1128

Air replacement:
Temperature: 45.0°C
Humidity: 20.1RH
Pressure: -1.00Pa
Sterilization time: 3hour
Air replace rounds: 1
round time: 6min(s)
2026-1-14 12:14

Air replacement:
Temperature: 45.5°C
Humidity: 20.1RH
Pressure: -1.00Pa
Sterilization time: 3hour
Air replace rounds: 2
round time: 6min(s)
2026-1-14 11:56

Air replacement:
Temperature: 45.9°C
Humidity: 20.1RH
Pressure: -1.00Pa
Sterilization time: 3hour
Air replace rounds: 3
round time: 6min(s)
2026-1-14 11:34

Air replacement:
Temperature: 46.0°C
Humidity: 20.1RH
Pressure: -1.00Pa
Sterilization time: 3hour
Air replace rounds: 4
round time: 6min(s)
2026-1-14 11:13

Air replacement:
Temperature: 45.2°C
Humidity: 20.1RH
Pressure: -1.00Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 6min(s)
2026-1-14 10:50

Sterilization:
Temperature: 45.2°C
Humidity: 20.1RH
Pressure: -0.00Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-14 10:47

Sterilization:
Temperature: 43.6°C
Humidity: 20.1RH
Pressure: -0.00Pa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-14 10:49

