

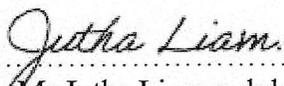
## Report Test

M/C: Henan Sanqiang Medical Equipment Liability Co. Ltd  
 Model: SQ-H220  
 Serial No: HNSQ2020050063

LOT No: 001131  
 Date: 20 January 2026

### Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
<b>Humidification</b>			
Humidity	50.0-90.0%	20.2 %	
<b>Sterilization</b>			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
<b>Gas Injection</b>			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
<b>Exposure</b>			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
<b>Exhaust</b>			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
<b>Storage</b>			
Storage Time	≥ 48 hour	≥ 48 hour	
<b>Chemical Test</b>			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
<b>Biological Test</b>			
EZ Test Biological Indicator	Green	Green	Pass

  
 (.....)  
 Mr. Jutha Liam mukda  
 Research and Development Engineer



			LOT	Color	BOX	TOTAL / PCS
			217525			
1	Z-INMLBSJ1-1 0-6 0-4-CR	STERILE SURGICAL STEEL INTERNALLY THREADED ROUND PRONG SET JEWELED MICRO LABRET	LOT NO 001131			4
2	Z-INMLBSJ1-1 0-8 0-4-CR	STERILE SURGICAL STEEL INTERNALLY THREADED ROUND PRONG SET JEWELED MICRO LABRET	LOT NO 001131			4
			217689			
1	Z-GPMLBJ-1 2-10 0-3-CR	STERILE GOLD PVD COATED SURGICAL STEEL JEWELED MICRO LABRET	LOT NO 001131			50
			217714			
2	Z-TINBN-3 0-19 0-6-HP	STERILE TITANIUM INTERNALLY THREADED CURVED BARBELL	LOT NO 001131			1
			217736			
1	Z-BNDJ-1 8-12 0-5-8-CR	STERILE SURGICAL STEEL DOUBLE JEWELED NAVEL BANANA	LOT NO 001131			100
2	Z-JNO-1 0-6 5-2 35-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELED CURVED NOSE STUD	LOT NO 001131			100
3	Z-MBSJ1-1 2-8 0-3-CR	STERILE SURGICAL STEEL ROUND PRONG SET JEWELED TRAGUS MICRO BARBELL	LOT NO 001131			50
4	Z-SES1-2 5-CR	STERILE SURGICAL STEEL ROUND PRONG SET JEWELED EAR STUDS PAIR	LOT NO 001131			50
			217742			
1	Z-TMBL-1 2-8 0-3-HP	STERILE TITANIUM MICRO BARBELL	LOT NO 001131			50
2	Z-TMCR-1 2-10 0-3-HP	STERILE TITANIUM MICRO CIRCULAR BARBELL	LOT NO 001131			50
			2563413			
4	THD-Z-BKNSJ148-1 2-CR	STERILE BLACK PVD COATED SURGICAL STEEL JEWELED MICRO ATTACHMENT FOR 1 2MM INTERNALLY THREADED PINS	LOT NO 001131			150
			2568458			
24	PT-Z-TLQNO-0-8 6-5-2-HP	STERILE TITANIUM LARGE CURVED BALL NOSE STUD	LOT NO 001131			150
38	PT-Z19-TBNDJ-1 8-10 0-5 8-CR	STERILE TITANIUM DOUBLE JEWELED NAVEL BANANA FOR PTIT ROC	LOT NO 001131		SET160	500
			2569293			
6	PT-Z-TBCH-1 2-12 0-4	STERILE TITANIUM BALL CLOSURE RING WITH HEJATITE BALL FOR PTIT ROC	LOT NO 001131			100
8	PT-Z-TBL-1 8-14 0-5	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 001131			200
10	PT-Z-TBL-1 6-22 0-5	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 001131			150
11	PT-Z-TBL-1 6-40 0-5	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 001131			100
15	PT-Z-TMBL-1 2-10 0-3-HP	STERILE TITANIUM MICRO BARBELL FOR PTIT ROC	LOT NO 001131			100
18	PT-Z-TMBN-1 2-10 0-3	STERILE TITANIUM MICRO CURVED BARBELL FOR PTIT ROC	LOT NO 001131			150
19	PT-Z-TMBH-1 2-8 0-3	STERILE TITANIUM MICRO CURVED BARBELL FOR PTIT ROC	LOT NO 001131			200
20	PT-Z-TMLB-1 2-10 0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 001131			500
21	PT-Z-TMLB-1 2-6 0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 001131			150
22	PT-Z-TMLB-1 2-7 0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 001131			100
24	PT-Z-TMLB-1 2-8 0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 001131			300
24	PT-Z-TMLB-1 2-8 0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 001131			250
33	PT-Z10-TMBN-1 2-10 0-3	STERILE TITANIUM MICRO CURVED BARBELL FOR PTIT ROC	LOT NO 001131		SET16	150
34	PT-Z10-TMBH-1 2-8 0-3	STERILE TITANIUM MICRO CURVED BARBELL FOR PTIT ROC	LOT NO 001131		SET16	150
			2567976			
4	SV-STERILIZE	SERVICE STERILE	LOT NO 001131			21
						3880

Sterilization:  
Temperature: 44.7°C  
Humidity: 20.1RH  
Pressures: -4.4kPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 15:22

Add ED gas  
Temperature: 44.5°C  
Humidity: 20.1RH  
Pressures: -4.0kPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 15:22

Evacuation  
Temperature: 44.2°C  
Humidity: 20.1RH  
Pressures: -0.0kPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 15:18

Keeping temp  
Temperature: 42.8°C  
Humidity: 20.1RH  
Pressures: -1.0kPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 14:18

Keep pressure qualified  
Temperature: 41.8°C  
Humidity: 20.1RH  
Pressures: -59.1kPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 14:15

Evacuation  
Temperature: 41.6°C  
Humidity: 20.1RH  
Pressures: -3.0kPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 14:14

Heating-up:  
Temperature: 23.5°C  
Humidity: 20.1RH  
Pressures: -3.0kPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 13:24

Setting:  
Temperature: 45.0°C  
Humidity: 55.0RH  
Pressures: -63.0kPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 13:24

20101126

END  
Lot 1131

Air replacement:  
Temperature: 45.4°C  
Humidity: 20.1RH  
Pressures: -1.0kPa  
Sterilization time: 0Hour  
Air replace rounds: 1  
round time: 0min(s)  
2026-1-20 19:50

Air replacement:  
Temperature: 45.4°C  
Humidity: 20.1RH  
Pressures: -1.0kPa  
Sterilization time: 0Hour  
Air replace rounds: 2  
round time: 0min(s)  
2026-1-20 19:29

Air replacement:  
Temperature: 45.5°C  
Humidity: 20.1RH  
Pressures: -1.0kPa  
Sterilization time: 0Hour  
Air replace rounds: 3  
round time: 0min(s)  
2026-1-20 19:9

Air replacement:  
Temperature: 45.5°C  
Humidity: 20.1RH  
Pressures: -1.0kPa  
Sterilization time: 0Hour  
Air replace rounds: 4  
round time: 0min(s)  
2026-1-20 18:49

Air replacement:  
Temperature: 45.1°C  
Humidity: 20.1RH  
Pressures: -1.0kPa  
Sterilization time: 0Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 16:23

Sterilization:  
Temperature: 44.3°C  
Humidity: 20.1RH  
Pressures: -1.4kPa  
Sterilization time: 0Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 16:23

Sterilization:  
Temperature: 44.1°C  
Humidity: 20.1RH  
Pressures: -5.1kPa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-1-20 15:24