

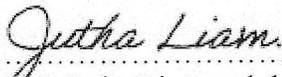
Report Test

M/C: Henan Sanqiang Medical Equipment Liability Co. Ltd
 Model: SQ-H220
 Serial No: HNSQ2020050063

LOT No: 001136
 Date: 28 January 2026

Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
Humidification			
Humidity	50.0-90.0%	20.2 %	
Sterilization			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
Gas Injection			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
Exposure			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
Exhaust			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
Storage			
Storage Time	≥ 48 hour	≥ 48 hour	
Chemical Test			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
Biological Test			
EZ Test Biological Indicator	Green	Green	Pass


 (.....)
 Mr. Jutha Liammukda
 Research and Development Engineer



		LOT	COBR	BOX	TOTAL / PCS
		217833			
3	Z-BNDJ-1.6-10.0-5/8-CR	STERILE SURGICAL STEEL DOUBLE JEWELLED NAVEL BANANA	LOT NO 001136		50
5	Z-XMBN-1.1-2-10-0-3-CL	STERILE BIOFLEX MICRO CURVED BARBELL WITH STEEL BALLS	LOT NO 001136		100
6	Z-XMBN-1.1-2-12-0-3-CL	STERILE BIOFLEX MICRO CURVED BARBELL WITH STEEL BALLS	LOT NO 001136		100
		217951			
1	Z-BC-1.0-10-0-4	STERILE SURGICAL STEEL BALL CLOSURE RING	LOT NO 001136		10
2	Z-BC-1.0-8-0-4	STERILE SURGICAL STEEL BALL CLOSURE RING	LOT NO 001136		5
3	Z-BC-1.2-10-0-4	STERILE SURGICAL STEEL BALL CLOSURE RING	LOT NO 001136		10
4	Z-BC-1.6-10-0-4	STERILE SURGICAL STEEL BALL CLOSURE RING	LOT NO 001136		10
5	Z-BC-2.0-10-0-4	STERILE SURGICAL STEEL BALL CLOSURE RING	LOT NO 001136		3
6	Z-BC-2.5-12-0-5	STERILE SURGICAL STEEL BALL CLOSURE RING	LOT NO 001136		3
7	Z-BL-1.6-10-0-4	STERILE SURGICAL STEEL BARBELL	LOT NO 001136		20
8	Z-BL-2.0-10-0-5	STERILE SURGICAL STEEL BARBELL	LOT NO 001136		3
9	Z-BL-2.5-12-0-5	STERILE SURGICAL STEEL BARBELL	LOT NO 001136		3
10	Z-BN-1.6-10-0-4	STERILE SURGICAL STEEL CURVED BARBELL	LOT NO 001136		10
11	Z-CB-1.6-8-0-3	STERILE SURGICAL STEEL CIRCULAR BARBELL	LOT NO 001136		10
12	Z-COBC-1.6-10-0	STERILE SURGICAL STEEL SEAMLESS RING	LOT NO 001136		3
13	Z-IDBLSC110-1.6-36.0-5-CR	STERILE SURGICAL STEEL INDUSTRIAL BARBELL WITH ADJUSTABLE SLIDING CHARM- FLOWER	LOT NO 001136		2
14	Z-LB-1.6-10-0-4	STERILE SURGICAL STEEL LABRET	LOT NO 001136		10
15	Z-MLB-1.2-10-0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 001136		5
16	Z-MLB-1.2-10-0-4	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 001136		10
17	Z-NO-1.0-6.5-2.5	STERILE SURGICAL STEEL 1.2MM THREADING CURVED NOSE STUD WITH BALL	LOT NO 001136		10
		217952			
7	Z-MBN-1.2-8-0-3	STERILE SURGICAL STEEL CURVED MICRO BARBELL	LOT NO 001136		10
8	Z-SES1-4-CR	STERILE SURGICAL STEEL ROUND PRONG SET JEWELLED EAR STUDS PAIR	LOT NO 001136		10
9	Z-SES1-4-JE	STERILE SURGICAL STEEL ROUND PRONG SET JEWELLED EAR STUDS PAIR	LOT NO 001136		5
		217954			
4	Z-SCES430-4	STERILE SURGICAL STEEL EAR STUDS PAIR - BALL 4MM	LOT NO 001136		1
		2570160			
2	PT-Z-GPTBL-1.6-20-0-5	STERILE GOLD PVD COATED TITANIUM BARBELL	LOT NO 001136		50
6	PT-Z-GPTMBL-1.2-10-0-3	STERILE GOLD PVD COATED TITANIUM MICRO BARBELL	LOT NO 001136		100
7	PT-Z-GPTMBN-1.2-8-0-3	STERILE GOLD PVD COATED TITANIUM CURVED MICRO BARBELL	LOT NO 001136		100
8	PT-Z-GPTMLB-1.2-8-0-3	STERILE GOLD PVD COATED TITANIUM MICRO LABRET	LOT NO 001136		400
9	PT-Z-GPTMLBJ-1.2-10-0-3-CR	STERILE GOLD PVD COATED TITANIUM PREMIUM CRYSTAL JEWELLED MICRO LABRET	LOT NO 001136		150
13	PT-Z-TBL-1.6-38-0-5	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 001136		50
14	PT-Z-TBNL-1.6-12-0-5/8-CN	SILHILE TITANIUM JEWELLED NAVEL BANANA FOR PTIII ROC	LOT NO 001136		50
15	PT-Z-TINBL-1.6-20-0-5-HP	STERILE TITANIUM INTERNALLY THREADED BARBELL FOR PTIT ROC	LOT NO 001136		50
16	PT-Z-TINMLB-1.2-8-0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO LABRET FOR PTIT ROC	LOT NO 001136		200
17	PT-Z-TMLB-1.2-12-0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 001136		100
18	PT-Z-TNO-0.8-6.5-2	STERILE TITANIUM CURVED BALL NOSE STUD FOR PTIT ROC	LOT NO 001136		150
22	PT-Z-TBL-1.6-16-0-5	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 001136	SET/25	250
25	PT-Z-TMLBJ-1.2-10-0-3-CR	STERILE TITANIUM MICRO LABRET WITH JEWELLED BALL FOR PTIT ROC	LOT NO 001136	SET/25	250
26	PT-Z-TMLBJ-1.2-8-0-3-CR	STERILE TITANIUM MICRO LABRET WITH JEWELLED BALL FOR PTIT ROC	LOT NO 001136	SET/20	200
		2567976			
8	SV-STERILIZE	SERVICE STERILE	LOT NO 001136		549
		2567982			
5	SV-STERILIZE	SERVICE STERILE	LOT NO 001136		4
		2567015			
4	SV-STERILIZE	SERVICE STERILE	LOT NO 001136		29
		2570404			
1	SV-STERILIZE	SERVICE STERILE	LOT NO 001136		5
					3050

Sterilization:
Temperature: 45.4°C
Humidity: 20.18H
Pressures: -4.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 15:23

Add EG gas
Temperature: 46.4°C
Humidity: 20.28H
Pressures: -4.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 15:29

Evacuation
Temperature: 44.9°C
Humidity: 20.18H
Pressures: -0.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 15:14

Keeping Temp
Temperature: 44.7°C
Humidity: 20.18H
Pressures: -1.6kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 14:14

Keep pressure qualified
Temperature: 43.1°C
Humidity: 20.18H
Pressures: -1.6kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 14:12

Evacuation
Temperature: 42.0°C
Humidity: 20.28H
Pressures: -3.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 14:0

Heating-up:
Temperature: 35.6°C
Humidity: 20.18H
Pressures: -3.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 13:38

Setting:
Temperature: 45.0°C
Humidity: 20.18H
Pressures: -6.0kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 13:36

28 Jan 26
END
Lot 1106

Air replacement:
Temperature: 45.1°C
Humidity: 20.28H
Pressures: -1.6kPa
Sterilization time: 3hour
Air replace rounds: 1
round time: 0min(s)
2026-1-28 19:50

Air replacement:
Temperature: 45.1°C
Humidity: 20.28H
Pressures: -1.6kPa
Sterilization time: 3hour
Air replace rounds: 2
round time: 0min(s)
2026-1-28 19:30

Air replacement:
Temperature: 46.3°C
Humidity: 20.18H
Pressures: -1.6kPa
Sterilization time: 3hour
Air replace rounds: 3
round time: 0min(s)
2026-1-28 19:10

Air replacement:
Temperature: 45.8°C
Humidity: 20.28H
Pressures: -1.6kPa
Sterilization time: 3hour
Air replace rounds: 4
round time: 0min(s)
2026-1-28 19:50

Air replacement:
Temperature: 45.1°C
Humidity: 20.18H
Pressures: -1.6kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 18:29

Sterilization:
Temperature: 45.0°C
Humidity: 20.18H
Pressures: -1.6kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 18:24

Sterilization:
Temperature: 45.7°C
Humidity: 20.18H
Pressures: -4.5kPa
Sterilization time: 3hour
Air replace rounds: 5
round time: 15min(s)
2026-1-28 15:25