

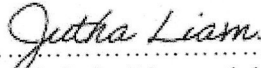
Report Test

M/C: Henan Sanqiang Medical Equipment Liability Co. Ltd
 Model: SQ-H220
 Serial No: HNSQ2020050063

LOT No: 001141
 Date: 04 February 2026

Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
Humidification			
Humidity	50.0-90.0%	20.2 %	
Sterilization			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
Gas Injection			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
Exposure			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
Exhaust			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
Storage			
Storage Time	≥ 48 hour	≥ 48 hour	
Chemical Test			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
Biological Test			
EZ Test Biological Indicator	Green	Green	Pass


 (.....)
 Mr. Jutha Liammukda
 Research and Development Engineer

			LOT	Color	BOX	TOTAL / PCS
			217523			
8	Z-GPTMBN-1.2-6.0-3	STERILE GOLD PVD COATED TITANIUM CURVED MICRO BARBELL	LOT NO 001141			3
1	Z-GBL-1.8-18.0-5	STERILE GOLD PVD COATED SURGICAL STEEL BARBELL	LOT NO 001141			50
2	Z-BNJ-1.6-10.0-6/8-CR	STERILE SURGICAL STEEL JEWELLED NAVAL BANANA	LOT NO 001141			100
1	Z-BNDJ-1.6-12.0-6/8-CR	STERILE SURGICAL STEEL DOUBLE JEWELLED NAVAL BANANA	LOT NO 001141			15
2	Z-MBN-1.2-10.0-3	STERILE SURGICAL STEEL CURVED MICRO BARBELL	LOT NO 001141			15
3	Z-MLB-1.2-8.0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 001141			15
4	Z-SCS430-3	STERILE SURGICAL STEEL EAR STUDS PAIR - BALL 3MM	LOT NO 001141			15
5	Z-BES1-2.5-CR	STERILE SURGICAL STEEL ROUND PRONG SET JEWELLED EAR STUDS PAIR	LOT NO 001141			20
6	Z-SCS1-3-CR	STERILE SURGICAL STEEL ROUND PRONG SET JEWELLED EAR STUDS PAIR	LOT NO 001141			10
7	Z-TNO-0.8-8.5-2-HP	STERILE TITANIUM CURVED BALL NOSE STUD				
			218124			
3	Z-BL-1.6-18.0-4	STERILE SURGICAL STEEL BARBELL	LOT NO 001141			26
4	Z-BNJ-1.6-12.0-5/8-CR	STERILE SURGICAL STEEL JEWELLED NAVAL BANANA	LOT NO 001141			20
			218125			
1	Z-BL-1.6-19.0-5	STERILE SURGICAL STEEL BARBELL	LOT NO 001141			10
2	Z-BL-1.6-40.0-4	STERILE SURGICAL STEEL BARBELL	LOT NO 001141			16
3	Z-MLB-1.2-10.0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 001141			15
4	Z-MLB-1.2-12.0-3	STERILE SURGICAL STEEL MICRO LABRET				
			218130			
1	Z-SES1-3-CR	STERILE SURGICAL STEEL ROUND PRONG SET JEWELLED EAR STUDS PAIR	LOT NO 001141			30
2	Z-LJNO-1.0-6.5-2.35-CR	STERILE SURGICAL STEEL PREMIUM CRYSTAL JEWELLED LARGE LEFT CURVE NOSE STUD	LOT NO 001141			40
3	Z-BNDJ-1.6-12.0-5/8-CR	STERILE SURGICAL STEEL DOUBLE JEWELLED NAVAL BANANA	LOT NO 001141			15
4	Z-MLB-1.2-8.0-3-CR	STERILE SURGICAL STEEL JEWELLED MICRO LABRET	LOT NO 001141			35
5	Z-MBN-1.2-11.2-3	STERILE SURGICAL STEEL CURVED MICRO BARBELL	LOT NO 001141			16
6	Z-BL-1.6-9.0-5	STERILE SURGICAL STEEL BARBELL	LOT NO 001141			10
7	Z-MCB-1.2-8.0-3	STERILE SURGICAL STEEL MICRO CIRCULAR BARBELL	LOT NO 001141			6
8	Z-BL-1.6-18.0-4	STERILE SURGICAL STEEL BARBELL	LOT NO 001141			10
9	Z-BL-1.6-14.0-4	STERILE SURGICAL STEEL BARBELL				
			218145			
9	Z-GBPCSH-0.8-8.0	GOLD PVD COATED SURGICAL STEEL HINGED SEGMENT RING PLAIN POLISHED	LOT NO 001141			30
10	Z-TBL-PIN-1.6-14.0-HP	STERILE TITANIUM BARBELL PIN	LOT NO 001141			50
11	Z-TBL-PIN-1.6-22.0-HP	STERILE TITANIUM BARBELL PIN	LOT NO 001141			50
12	Z-TBN-PIN-1.6-10.0-HP	STERILE TITANIUM CURVED BARBELL PIN	LOT NO 001141			30
13	Z-TBCSH-0.8-8.0-HP	STERILE TITANIUM HINGED SEGMENT RING	LOT NO 001141			50
14	Z-TBCSH-1.0-7.0-HP	STERILE TITANIUM HINGED SEGMENT RING	LOT NO 001141			50
15	Z-TBCSH-1.0-9.0-HP	STERILE TITANIUM HINGED SEGMENT RING				
			218146			
5	Z-TBL-PIN-1.6-18.0-HP	STERILE TITANIUM BARBELL PIN	LOT NO 001141			50
6	Z-TBL-PIN-1.6-18.0-HP	STERILE TITANIUM BARBELL PIN	LOT NO 001141			50
			218149			
2	Z-BNDJ-1.6-10.0-4/8-CR	STERILE SURGICAL STEEL DOUBLE PREMIUM CRYSTALS JEWELLED NAVAL BANANA	LOT NO 001141			5
3	Z-BNDJ-1.6-12.0-4/8-CR	STERILE SURGICAL STEEL DOUBLE PREMIUM CRYSTALS JEWELLED NAVAL BANANA	LOT NO 001141			10
4	Z-MBNF4-1.2-8.0-3	STERILE SURGICAL STEEL TRIPLE STARS FANCY CURVED MICRO BARBELL	LOT NO 001141			3
5	Z-TINBNDJ-1.6-10.0-5/8-HP-CR	STERILE TITANIUM INTERNALLY THREADED PREMIUM CRYSTAL JEWELLED NAVAL BANANA	LOT NO 001141			2
6	Z-TLJNO-S-0.8-19.0-1.8-HP-CR	STERILE TITANIUM PREMIUM CRYSTAL JEWELLED STRAIGHT NOSE STUD 18MM	LOT NO 001141			5
			218151			
1	Z-TNO-S-0.8-15.0-2-HP	STERILE TITANIUM STRAIGHT BALL NOSE STUD 16MM	LOT NO 001141			2
			218157			
5	Z-QER007-0.8-8.0	STERILE SURGICAL STEEL OPEN CLAW EARRING	LOT NO 001141			12
			2570335			
9	THD-Z-BKTINMC-1.2-3	STERILE BLACK PVD COATED TITANIUM MICRO CONE FOR 1.2MM INTERNALLY THREADED PINS	LOT NO 001141			100
10	THD-Z-BKTINPD2-1.6-4	BLACK PVD COATED TITANIUM GLITTERLINE DISC FOR 1.6MM INTERNALLY THREADED PINS	LOT NO 001141			100
37	THD-Z-MLB-1.2-10.0-3	STERILIZED SURGICAL STEEL MICRO LABRET FOR THENDRO	LOT NO 001141			1200
40	THD-Z-SBNDJ-1.6-10.0-4/8-CR	STERILE SURGICAL STEEL DOUBLE JEWELLED MINI NAVAL BANANA FOR THENDRO	LOT NO 001141			300
42	THD-Z-TINBL-1.6-8.0-4	STERILE TITANIUM INTERNALLY THREADED CURVED BARBELL	LOT NO 001141			500
43	THD-Z-TINBN-1.6-10.0-4	STERILE TITANIUM INTERNALLY THREADED CURVED BARBELL	LOT NO 001141			300
54	THD-Z-TINMLB-1.2-8.0-3	STERILE TITANIUM INTERNALLY THREADED MICRO LABRET	LOT NO 001141			2500
54	THD-Z-TINMLB-1.2-8.0-3	STERILE ANODISED TITANIUM INTERNALLY THREADED JEWELLED MICRO LABRET	LOT NO 001141			200
55	THD-Z-TINPD2-1.6-4-HP	TITANIUM GLITTERLINE DISC FOR 1.6MM INTERNALLY THREADED PINS	LOT NO 001141			100
55	THD-Z-TINPD2-1.6-4-HP	TITANIUM GLITTERLINE DISC FOR 1.6MM INTERNALLY THREADED PINS	LOT NO 001141			400
58	THD-Z-TMLB-1.2-10.0-3	STERILE TITANIUM MICRO CURVED BARBELL FOR THENDRO	LOT NO 001141			600
58	THD-Z-TMLB-1.2-8.0-3	STERILE TITANIUM MICRO LABRET FOR THENDRO				
			2570836			
1	Z-TINMLB-PIN-1.2-8.0-3B-HP	STERILE TITANIUM INTERNALLY THREADED MICRO LABRET PIN FOR 0.9MM THREAD	LOT NO 001141			10
			2570862			
45	THD-Z-TINJD-1.6-4-RN-CR	STERILE GOLD ANODISED TITANIUM JEWELLED DISC FOR 1.6MM INTERNALLY THREADED PINS	LOT NO 001141			100
			2570867			
5	PT-Z-GPTNO-0.8-6.5-2	STERILE GOLD PVD COATED TITANIUM CURVED BALL NOSE STUD	LOT NO 001141			150
6	PT-Z-TBCSH-0.8-8.0-HP	STERILE TITANIUM HINGED SEGMENT RING FOR PTIT ROC	LOT NO 001141			60
7	PT-Z-TBN-1.6-10.0-4	STERILE TITANIUM TITANIUM CURVED BARBELL FOR PTIT ROC	LOT NO 001141			30
11	PT-Z-TJNO-0.8-6.5-2.35-CR	STERILE TITANIUM JEWELLED CURVED NOSE STUD FOR PTIT ROC	LOT NO 001141		SET/15	150
18	PT-Z10-TMLBJ-1.2-8.0-2.5-CR	STERILE TITANIUM MICRO LABRET WITH JEWELLED BALL FOR PTIT ROC	LOT NO 001141			
			2567982			
1	SV-STERILIZE	SERVICE STERILE	LOT NO 001141			4
						8205

Sterilization:
Temperature: 44.8°C
Humidity: 20.1RH
Pressure: -4.1kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 13:54

Add ED gas
Temperature: 44.7°C
Humidity: 20.1RH
Pressure: -4.0kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 13:54

Evacuation
Temperature: 44.1°C
Humidity: 20.1RH
Pressure: -9.0kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 13:38

Keeping Temp
Temperature: 42.3°C
Humidity: 20.1RH
Pressure: -1.6kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 12:38

Keep pressure qualified
Temperature: 41.2°C
Humidity: 20.1RH
Pressure: -59.3kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 12:36

Evacuation
Temperature: 41.4°C
Humidity: 20.1RH
Pressure: -8.0kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 12:25

Heating up
Temperature: 21.8°C
Humidity: 20.1RH
Pressure: -3.0kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 11:41

Setting
Temperature: 45.0°C
Humidity: 55.0RH
Pressure: -60.0kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 11:40

04/02/20
END
Lot 1141

Air replacement:
Temperature: 45.0°C
Humidity: 20.1RH
Pressure: -1.6kPa
Sterilization time: 3Hour
Air replace rounds: 1
round time: 0min(s)
2020-2-4 16:21

Air replacement:
Temperature: 45.3°C
Humidity: 20.2RH
Pressure: -1.6kPa
Sterilization time: 3Hour
Air replace rounds: 2
round time: 0min(s)
2020-2-4 16:11

Air replacement:
Temperature: 45.6°C
Humidity: 20.1RH
Pressure: -1.6kPa
Sterilization time: 3Hour
Air replace rounds: 3
round time: 0min(s)
2020-2-4 17:41

Air replacement:
Temperature: 45.7°C
Humidity: 20.1RH
Pressure: -1.6kPa
Sterilization time: 3Hour
Air replace rounds: 4
round time: 0min(s)
2020-2-4 17:20

Air replacement:
Temperature: 45.1°C
Humidity: 20.1RH
Pressure: -1.6kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 17:00

Sterilization:
Temperature: 44.6°C
Humidity: 20.2RH
Pressure: -2.4kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 16:54

Sterilization:
Temperature: 45.0°C
Humidity: 20.1RH
Pressure: -4.4kPa
Sterilization time: 3Hour
Air replace rounds: 5
round time: 15min(s)
2020-2-4 13:56

