

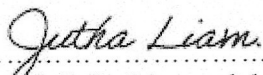
Report Test

M/C: Henan Sanqiang Medical Equipment Liability Co. Ltd
 Model: SQ-H220
 Serial No: HNSQ2020050063

LOT No: 001163
 Date: 10 March 2026

Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
Humidification			
Humidity	50.0-90.0%	20.2 %	
Sterilization			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
Gas Injection			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
Exposure			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
Exhaust			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
Storage			
Storage Time	≥ 48 hour	≥ 48 hour	
Chemical Test			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
Biological Test			
EZ Test Biological Indicator	Green	Green	Pass


 (.....)
 Mr. Jutha Liammukda
 Research and Development Engineer



			LOT	Color	BOX	TOTAL / PCS
			219036			
1	Z-TBL-1.6-38.0-4-HP	STERILE TITANIUM BARBELL	LOT NO 001163			5
2	Z-TBNDJ-1.6-12.0-5/8-HP-CR	STERILE TITANIUM DOUBLE JEWELLED NAVEL BANANA	LOT NO 001163			5
3	Z-TPFINJD-1.0-8.0-3-HP-CR	STERILE TITANIUM JEWELLED THREADLESS DISC LABRET	LOT NO 001163			8
4	Z-TPFINJD-1.2-10.0-3-HP-CR	STERILE TITANIUM JEWELLED THREADLESS DISC LABRET	LOT NO 001163			6
			219041			
1	Z-TINBL-1.6-22.0-5-HP	STERILE TITANIUM INTERNALLY THREADED BARBELL	LOT NO 001163			15
2	Z-TINDAYSJD-1.6-2.5-4-HP-CR	STERILE TITANIUM INTERNALLY THREADED DERMAL ANCHOR BIG HOLE WITH JEWELLED DISC	LOT NO 001163			6
3	Z-TINMBN-1.2-10.0-3-HP	STERILE TITANIUM INTERNALLY THREADED CURVED MICRO BARBELL	LOT NO 001163			20
4	Z-TINMCB-1.2-10.0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO CIRCULAR BARBELL	LOT NO 001163			20
5	Z-TINMLB-1.2-10.0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO LABRET	LOT NO 001163			20
6	ZDPU-1.5	STERILE DERMAL PUNCH	LOT NO 001163			20
7	ZHM	STERILE TWEEZERS 5.5 INCHES LONG	LOT NO 001163			5
8	ZUPLST-10.0-15.0	STERILE POLYMER PIERCING CORK STOPPER	LOT NO 001163			20
			219085			
1	Z-TINMBL-1.2-12.0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO BARBELL	LOT NO 001163			25
			219086			
1	Z-TINMBL-1.2-8.0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO BARBELL	LOT NO 001163			25
2	Z-TINMBN-1.2-12.0-3-HP	STERILE TITANIUM INTERNALLY THREADED CURVED MICRO BARBELL	LOT NO 001163			25
3	Z-TINMLBJ-1.2-8.0-2-HP-CR	STERILE TITANIUM INTERNALLY THREADED JEWELLED MICRO LABRET	LOT NO 001163			50
			2569779			
13	TDI-Z-TBNJ-OP-1.6-12.0-5/8-CR	STERILE TITANIUM JEWELLED NAVEL BANANA FOR TDI	LOT NO 001163			100
			2574095			
35	PT-Z-XBL-1-1.6-18.0-5-CL	STERILE BIOFLEX BARBELL WITH STEEL BALLS FOR PTIT ROC	LOT NO 001163			25
36	PT-Z-XMBN-1-1.2-10.0-3-CL	STERILE BIOFLEX MICRO CURVED BARBELL WITH STEEL BALLS FOR PTIT ROC	LOT NO 001163			20
38	PT-Z-XMLB-1-1.2-10.0-3-4B-CL	STERILE BIOFLEX MICRO LABRET WITH STEEL BALL FOR PTIT ROC	LOT NO 001163			50
			2574032			
12	THD-Z-GPMCB-1.2-10.0-3	STERILE GOLD PVD COATED SURGICAL STEEL MICRO CIRCULAR BARBEL FOR THENDRO	LOT NO 001163			200
19	THD-Z-MBN-1.2-10.0-3	STERILE SURGICAL STEEL MICRO CURVED BARBELL FOR THENDRO	LOT NO 001163			350
21	THD-Z-MCB-1.2-8.0-3	STERILE SURGICAL STEEL MICRO CIRCULAR BARBELL FOR THENDRO	LOT NO 001163			200
22	THD-Z-MLB-1.2-10.0-3	STERILIZED SURGICAL STEEL MICRO LABRET FOR THENDRO	LOT NO 001163			700
29	THD-Z-TINBL-1.6-20.0-5	STERILE TITANIUM INTERNALLY THREADED BARBELL	LOT NO 001163			250
32	THD-Z-TINLB-1.6-8.0-4	STERILE TITANIUM INTERNALLY THREADED LABRET	LOT NO 001163			200
33	THD-Z-TINMBL-1.2-8.0-3	STERILE TITANIUM INTERNALLY THREADED MICRO BARBELL	LOT NO 001163			1000
41	THD-Z-TIMBN-1.2-8.0-3	STERILE TITANIUM MICRO CURVED BARBELL FOR THENDRO	LOT NO 001163			250
42	THD-Z-TMLB-1.2-10.0-3	STERILE TITANIUM MICRO LABRET FOR THENDRO	LOT NO 001163			800
			2574630			
1	SV-STERILIZE	SERVICE STERILE	LOT NO 001163			2
			2570610			
1	SV-STERILIZE	SERVICE STERILE	LOT NO 001163			253
						4675

Sterilization:
 Temperature: 44.7°C
 Humidity: 20.28%
 Pressures: -1.6kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 13:39

Add E9 gas
 Temperature: 44.5°C
 Humidity: 20.18%
 Pressures: -4.0kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 13:38

Evacuation
 Temperature: 44.2°C
 Humidity: 20.18%
 Pressures: -0.0kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 13:34

Keeping Temp
 Temperature: 42.6°C
 Humidity: 20.18%
 Pressures: -1.6kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 12:54

Keep pressure qualified
 Temperature: 41.4°C
 Humidity: 20.18%
 Pressures: -53.0kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 12:32

Evacuation
 Temperature: 41.4°C
 Humidity: 20.28%
 Pressures: -3.0kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 12:29

Heating-up:
 Temperature: 21.5°C
 Humidity: 20.18%
 Pressures: -3.0kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 11:38

Setting:
 Temperature: 45.0°C
 Humidity: 25.08%
 Pressures: -60.0kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 11:35

No 10/3/26
 Lot 1167

Air replacement:
 Temperature: 45.9°C
 Humidity: 20.28%
 Pressures: -1.6kPa
 Sterilization time: 3hour
 Air replace rounds: 1
 round time: 0min(s)
 2026-3-10 18:17

Air replacement:
 Temperature: 45.5°C
 Humidity: 20.28%
 Pressures: -1.0kPa
 Sterilization time: 3hour
 Air replace rounds: 2
 round time: 0min(s)
 2026-3-10 17:47

Air replacement:
 Temperature: 45.0°C
 Humidity: 20.18%
 Pressures: -1.6kPa
 Sterilization time: 3hour
 Air replace rounds: 3
 round time: 0min(s)
 2026-3-10 17:26

Air replacement:
 Temperature: 45.6°C
 Humidity: 20.18%
 Pressures: -1.6kPa
 Sterilization time: 3hour
 Air replace rounds: 4
 round time: 0min(s)
 2026-3-10 17:16

Air replacement:
 Temperature: 45.6°C
 Humidity: 20.18%
 Pressures: -1.0kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 16:45

Sterilization:
 Temperature: 44.5°C
 Humidity: 20.18%
 Pressures: -0.0kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 16:39

Sterilization:
 Temperature: 44.1°C
 Humidity: 20.18%
 Pressures: -0.0kPa
 Sterilization time: 3hour
 Air replace rounds: 5
 round time: 15min(s)
 2026-3-10 13:41