

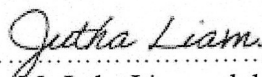
## Report Test

M/C: Henan Sanqiang Medical Equipment Liability Co. Ltd  
 Model: SQ-H220  
 Serial No: HNSQ2020050063

LOT No: 001173  
 Date: 25 March 2026

### Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
<b>Humidification</b>			
Humidity	50.0-90.0%	20.2 %	
<b>Sterilization</b>			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
<b>Gas Injection</b>			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
<b>Exposure</b>			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
<b>Exhaust</b>			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
<b>Storage</b>			
Storage Time	≥ 48 hour	≥ 48 hour	
<b>Chemical Test</b>			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
<b>Biological Test</b>			
EZ Test Biological Indicator	Green	Green	Pass

  
 (.....)  
 Mr. Jutha Liammukda  
 Research and Development Engineer



			LOT	Color	BOX	TOTAL / PCS
			219489			
1	Z-BL-1 6-19 0-5	STERILE SURGICAL STEEL BARBELL	LOT NO 1173			100
2	Z-MBL-1 2-10 0-3	STERILE SURGICAL STEEL MICRO BARBELL	LOT NO 1173			200
3	Z-MLB-1 2-8 0-3	STERILE SURGICAL STEEL MICRO LABRET	LOT NO 1173			150
			219531			
1	Z-TINMCB-1 2-10 0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO CIRCULAR BARBELL	LOT NO 1173			12
2	Z-TINMBN-1 2-12 0-3-HP	STERILE TITANIUM INTERNALLY THREADED CURVED MICRO BARBELL	LOT NO 1173			12
4	Z-TINBNGU-1 8-12 0-5-8-HP-CR	STERILE TITANIUM INTERNALLY THREADED PREMIUM CRYSTAL JEWELLED NAVEL BANANA	LOT NO 1173			12
			219586			
2	Z-GPSES1-2 5-CR	STERILE GOLD PVD COATED SURGICAL STEEL ROUND PRONG SET JEWELLED EAR STUDS PAIR	LOT NO 1173			5
			219575			
1	Z-MCB-1 2-10 0-3	STERILE SURGICAL STEEL MICRO CIRCULAR BARBELL	LOT NO 1173			100
2	Z-TINBU-1 6-12 0-5-8-HP-CR	STERILE TITANIUM DOUBLE JEWELLED NAVEL BANANA	LOT NO 1173			50
3	Z-TINCO-1 0-6 5-2 35-HP-CR	STERILE TITANIUM PREMIUM CRYSTAL JEWELLED CURVED NOSE STUD	LOT NO 1173			41
4	Z-TMLBJ-1 2-8 0-3-HP-CR	TITANIUM GRADE 23 LABRET FLAT BACK SINGLE GEM BALL END	LOT NO 1173			50
			2568458			
5	PT-Z-GPTSNJ-1 8-12 0-5-8-CR	STERILE GOLD PVD COATED TITANIUM PREMIUM CRYSTAL JEWELLED NAVEL BANANA	LOT NO 1173			150
			2570393			
39	THD-Z-RGNSC288-1 2	STERILE ROSE GOLD PVD COATED SURGICAL STEEL MICRO ATTACHMENT FOR 1.2MM INTERNALLY THREADED PINS - SWORD	LOT NO 1173			100
			2574095			
5	PT-Z-GPTMLBJ-1 2-8 0-3-CR	STERILE GOLD PVD COATED TITANIUM PREMIUM CRYSTAL JEWELLED MICRO LABRET	LOT NO 1173			250
			2575610			
2	PT-Z-GPTBL-1 8-18 0-5	STERILE GOLD PVD COATED TITANIUM BARBELL	LOT NO 1173			100
			2575841			
26	PT-Z-XBINJ-1 1 6-12 0-5-8-CR	STERILE BIOFLEX JEWELLED NAVEL BANANA WITH SURGICAL STEEL BALL FOR PTIT ROC	LOT NO 1173			50
32	PT-Z10-TMLBJ-1 2-8 0-2 5-CR	STERILE TITANIUM MICRO LABRET WITH JEWELLED BALL FOR PTIT ROC	LOT NO 1173		SET115	150
			2576113			
1	ZPRTA-6 0-70 0	STERILE DISPOSABLE MULTI ANGLED RECEIVING TUBE 8 0MM	LOT NO 1173			1558
			2576247			
1	TDI-Z-TBL-1 6-14 0-5	STERILE TITANIUM BARBELL FOR TDI	LOT NO 1173			300
2	TDI-Z-TBL-1 8-16 0-5	STERILE TITANIUM BARBELL FOR TDI	LOT NO 1173			300
5	TDI-Z-TBL-1 6-22 0-5	STERILE TITANIUM BARBELL FOR TDI	LOT NO 1173			300
			2576395			
1	Z-TINMLBJ-1 2-10 0-2-HP-CR	STERILE TITANIUM INTERNALLY THREADED JEWELLED MICRO LABRET	LOT NO 1173			15
			2575000			
2	SV-STERILIZE	SERVICE STERILE	LOT NO 1173			3
			2574865			
1	SV-STERILIZE	SERVICE STERILE	LOT NO 1173			2
			2574997			
1	SV-STERILIZE	SERVICE STERILE	LOT NO 1173			10
			2576412			
1	SV-STERILIZE	SERVICE STERILE	LOT NO 1173			25
						4045

Sterilization:  
Temperature: 43.8°C  
Humidity: 20.18%  
Pressure: -0.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-3-25 13:36

Add GA gas  
Temperature: 43.6°C  
Humidity: 20.18%  
Pressure: -4.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-3-25 13:36

Evacuation  
Temperature: 44.2°C  
Humidity: 20.18%  
Pressure: -0.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-3-25 13:31

Keeping Temp  
Temperature: 42.8°C  
Humidity: 20.18%  
Pressure: -1.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-3-25 12:31

Keep pressure qualified  
Temperature: 41.6°C  
Humidity: 20.18%  
Pressure: -0.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-3-25 12:29

Evacuation  
Temperature: 41.5°C  
Humidity: 20.28%  
Pressure: -3.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-3-25 12:11

Heating-up  
Temperature: 32.3°C  
Humidity: 20.18%  
Pressure: -3.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-3-25 11:37

Setting:  
Temperature: 45.0°C  
Humidity: 20.18%  
Pressure: -0.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-3-25 11:30

25/3/26  
END Lot 1173

Air replacement:  
Temperature: 45.2°C  
Humidity: 20.28%  
Pressure: -1.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 1  
round time: 0min(s)  
2026-3-25 18:39

Air replacement:  
Temperature: 47.0°C  
Humidity: 20.18%  
Pressure: -1.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 2  
round time: 0min(s)  
2026-3-25 18:14

Air replacement:  
Temperature: 45.6°C  
Humidity: 20.18%  
Pressure: -2.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 3  
round time: 0min(s)  
2026-3-25 17:37

Air replacement:  
Temperature: 46.4°C  
Humidity: 20.18%  
Pressure: -1.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 4  
round time: 0min(s)  
2026-3-25 17:14

Air replacement:  
Temperature: 45.8°C  
Humidity: 20.18%  
Pressure: -1.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 0min(s)  
2026-3-25 16:48

Sterilization:  
Temperature: 45.0°C  
Humidity: 20.28%  
Pressure: -0.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-3-25 16:38

Sterilization:  
Temperature: 45.7°C  
Humidity: 20.18%  
Pressure: -0.06Pa  
Sterilization time: 3Hour  
Air replace rounds: 5  
round time: 15min(s)  
2026-3-25 13:38