


## Report Test


M/C: Henan Sanqiang Medical Equipment Liability Co. Ltd  
 Model: SQ-H220  
 Serial No: HNSQ2020050063

LOT No: 001183  
 Date: 21 April 2026

### Determination of Sterilization Operation Specification

Item Detail	Specification	Setting	Result
Chamber Setting Temperature	40.0-55.5 °C	45.0 °C	
Vacuum pressure	40-60 kPa	60 kPa	
Vacuum Time	50-60 min	60 min	
<b>Humidification</b>			
Humidity	50.0-90.0%	20.2 %	
<b>Sterilization</b>			
Sterilization Temperature	40.0-55.5 °C	43.5 °C	
<b>Gas Injection</b>			
Final pressure	0-10 kPa	4 kPa	
Final temperature of vaporizer	40.0-55.5 °C	44.8 °C	
<b>Exposure</b>			
Temperature at start point of Exposure time period	40.0-55.5 °C	43.3 °C	
Temperature at finish point of Exposure time period	40.0-55.5 °C	46.0 °C	
Exposure Time	120-240 min	180 min	
<b>Exhaust</b>			
Exhaust Time	60-80 min	75 min	
First Vacuum pressure	0-2 kPa	1.5 kPa	
Second Vacuum pressure	0-2 kPa	1.5 kPa	
Third Vacuum pressure	0-2 kPa	1.5 kPa	
Fourth Vacuum pressure	0-2 kPa	1.5 kPa	
Fifth Vacuum pressure	0-2 kPa	1.5 kPa	
<b>Storage</b>			
Storage Time	≥ 48 hour	≥ 48 hour	
<b>Chemical Test</b>			
EO sterilization temperature 40-55°C, Time Sterilization by EO 3hours	Blue	Red	Pass
<b>Biological Test</b>			
EZ Test Biological Indicator	Green	Green	Pass

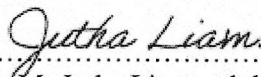  
 (.....)  
 Mr. Jutha Liammukda  
 Research and Development Engineer

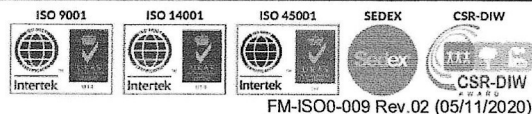

FM-ISO0-009 Rev.02 (05/11/2020)

			LOT	Cabr	BOX	TOTAL / PCS
		220018				
1	Z-BNDJ-1.6-14.0-5/8-CR	STERILE SURGICAL STEEL DOUBLE PREMIUM CRYSTALS JEWELED NAVEL BANANA	LOT NO 001183			100
2	Z-MBNDJ-1.2-8.0-3-CR	STERILE SURGICAL STEEL DOUBLE JEWELED CURVED MICRO BARBELL	LOT NO 001183			125
		220043				
1	Z-TINMLBJ-1.0-6.0-2-HP-CR	STERILE TITANIUM INTERNALLY THREADED JEWELED MICRO LABRET	LOT NO 001183			20
2	Z-TINMLBJ-1.2-7.0-2-HP-CR	STERILE TITANIUM INTERNALLY THREADED JEWELED MICRO LABRET	LOT NO 001183			20
3	Z-TINMLBJ-1.2-8.0-3-HP-CR	STERILE TITANIUM INTERNALLY THREADED JEWELED MICRO LABRET	LOT NO 001183			20
		220127				
1	Z-TESS3-3-CR	STERILE PAIR OF TITANIUM SQUARE PRONG SET JEWELED EAR STUDS	LOT NO 001183			10
		220142				
2	ZDPU-1.5	STERILE DERMAL PUNCH	LOT NO 001183			10
		2576385				
1	TDI-Z-TINMBN-1.2-10.0-3	STERILE TITANIUM INTERNALLY THREADED CURVED MICRO BARBELL FOR TDI	LOT NO 001183			100
		2578537				
1	TDI-Z-TMBL-1.2-8.0-3-HP	STERILE TITANIUM MICRO BARBELL FOR TDI	LOT NO 001183			100
		2578314				
4	PT-Z-GPTMCB-1.2-10.0-3	STERILE GOLD PVD COATED TITANIUM MICRO CIRCULAR BARBELL	LOT NO 001183			100
5	PT-Z-GPTMCB-1.2-8.0-3	STERILE GOLD PVD COATED TITANIUM MICRO CIRCULAR BARBELL	LOT NO 001183			100
6	PT-Z-GPTMLB-1.2-10.0-3	STERILE GOLD PVD COATED TITANIUM MICRO LABRET	LOT NO 001183			300
9	PT-Z-GPTNO-0.8-6.5-2	STERILE GOLD PVD COATED TITANIUM CURVED BALL NOSE STUD	LOT NO 001183			150
11	PT-Z-TBL-1.6-12.0-4	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 001183			50
12	PT-Z-TBL-1.6-40.0-5	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 001183			50
13	PT-Z-TBN-1.6-10.0-4	STERILE TITANIUM TITANIUM CURVED BARBELL FOR PTIT ROC	LOT NO 001183			50
15	PT-Z-TINMLB-1.2-8.0-3-HP	STERILE TITANIUM INTERNALLY THREADED MICRO LABRET FOR PTIT ROC	LOT NO 001183			200
16	PT-Z-TMBL-1.2-10.0-3-HP	STERILE TITANIUM MICRO BARBELL FOR PTIT ROC	LOT NO 001183			100
18	PT-Z-TMBN-1.2-14.0-3	STERILE TITANIUM MICRO CURVED BARBELL FOR PTIT ROC	LOT NO 001183			50
19	PT-Z-TMLB-1.2-12.0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 001183			100
25	PT-Z10-TBL-1.6-14.0-5	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 001183		SET/20	200
26	PT-Z10-TBL-1.6-18.0-5	STERILIZED TITANIUM BARBELL FOR PTIT ROC	LOT NO 001183		SET/15	150
29	PT-Z10-TMBN-1.2-10.0-3	STERILE TITANIUM MICRO CURVED BARBELL FOR PTIT ROC	LOT NO 001183		SET/20	200
30	PT-Z10-TMLB-1.2-8.0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 001183		SET/100	1000
31	PT-Z10-TMLB-1.2-9.0-3	STERILE TITANIUM MICRO LABRET FOR PTIT ROC	LOT NO 001183		SET/50	500
32	PT-Z10-TMLBJ-1.2-10.0-2.5-CR	STERILE TITANIUM MICRO LABRET WITH JEWELED BALL FOR PTIT ROC	LOT NO 001183		SET/15	150
		2576242				
2	SV-STERILIZE	SERVICE STERILE	LOT NO 001183			15
		2577508				
5	SV-STERILIZE	SERVICE STERILE	LOT NO 001183			1
		2578284				
1	SV-STERILIZE	SERVICE STERILE	LOT NO 001183			3
		2578649				
1	SV-STERILIZE	SERVICE STERILE	LOT NO 001183			24
		2578772				
2	SV-STERILIZE	SERVICE STERILE	LOT NO 001183			40
						4038

Sterilization:  
Temperature: 44.5°C  
Humidity: 20.19%  
Pressures: -2.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 15min(s)  
2026-4-21 13:24

Add ED gas  
Temperature: 44.5°C  
Humidity: 20.19%  
Pressures: -4.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 15min(s)  
2026-4-21 13:24

Evacuation  
Temperature: 44.7°C  
Humidity: 20.29%  
Pressures: -0.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 15min(s)  
2026-4-21 13:19

Keeping Temp  
Temperature: 43.8°C  
Humidity: 20.19%  
Pressures: -1.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 15min(s)  
2026-4-21 12:19

Feed pressure qualified  
Temperature: 41.7°C  
Humidity: 20.19%  
Pressures: -0.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 15min(s)  
2026-4-21 12:17

Evacuation  
Temperature: 41.5°C  
Humidity: 20.29%  
Pressures: -0.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 15min(s)  
2026-4-21 12:15

Heating on  
Temperature: 20.9°C  
Humidity: 20.19%  
Pressures: -3.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 15min(s)  
2026-4-21 11:26

Settling  
Temperature: 45.0°C  
Humidity: 20.69%  
Pressures: -0.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 15min(s)  
2026-4-21 11:26

21104126  
Lab 1183  
END

Air replacement:  
Temperature: 45.8°C  
Humidity: 20.29%  
Pressures: -1.0kPa  
Sterilization time: 30min  
Air replace rounds: 1  
round time: 0min(s)  
2026-4-21 17:54

Air replacement:  
Temperature: 45.2°C  
Humidity: 20.29%  
Pressures: -1.0kPa  
Sterilization time: 30min  
Air replace rounds: 2  
round time: 0min(s)  
2026-4-21 17:34

Air replacement:  
Temperature: 45.3°C  
Humidity: 20.19%  
Pressures: -1.0kPa  
Sterilization time: 30min  
Air replace rounds: 3  
round time: 0min(s)  
2026-4-21 17:19

Air replacement:  
Temperature: 45.5°C  
Humidity: 20.29%  
Pressures: -1.0kPa  
Sterilization time: 30min  
Air replace rounds: 4  
round time: 0min(s)  
2026-4-21 16:52

Air replacement:  
Temperature: 45.9°C  
Humidity: 20.29%  
Pressures: -1.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 0min(s)  
2026-4-21 16:31

Sterilization:  
Temperature: 44.4°C  
Humidity: 20.29%  
Pressures: -0.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 15min(s)  
2026-4-21 16:25

Sterilization:  
Temperature: 44.2°C  
Humidity: 20.29%  
Pressures: -0.0kPa  
Sterilization time: 30min  
Air replace rounds: 5  
round time: 15min(s)  
2026-4-21 15:26